



For a good experience, please read and keep this manual carefully before use.

1. Press the power key for a long time to boot for 2 seconds, and you can't use it on when charging.
2. Be sure to remove the jewelry (necklaces, earrings, etc.) (Prevent jewelry from being caught by the Roller Massager).
3. Please tie up the hair (the hair may be caught in it, and the hair around the corner and the short hair should be fixed and used. When the roller Massager is wrapped around the hair, Please removing the hair before turning it off
4. When using only the Roller Massager (OFF), the roller massager can be directly affixed to the skin Use.
5. When the Roller Massager is used with EMS function (ON), please use water or makeup water, gel and other moisturizing skin before using it.
6. At the same time, two rollers contact the skin to form a loop in order to feel the stimulation of EMS.
7. This product uses micro-current ion technology, if there is some electric flu is a normal phenomenon, please rest assured to use.
8. When using, do not stay in the same area for more than three seconds to avoid skin discomfort.


Roller Massager is a beauty instrument which enhances cell activity, desalinizes wrinkles, elevates facial outline, and integrates eccentric roller + temperature and heat + micro current.

- * 3D texture roller: eccentric shaft roller to achieve manual massage push, pull, kneading techniques.
- * Mild: warm massage promotes absorption, relaxes the texture of the skin, and makes the skin fuller.
- * Microcurrent: use voltage stimulation to tighten human skin, increase its own elasticity, and enhance the metabolism of skin cells.

1. Please pay attention to adjust the angle between the roller and the skin during use to avoid pinching the skin. Do not use it when the skin is tightly clamped.
2. Do not rub hard on the skin, especially in the neck or arm, please pay attention to the intensity of nursing when using.
3. When the roller is used with EMS function (turn on), please use water or makeup water, gel and other moisturizing skin before using it.
4. In order to feel the stimulation of EMS, the two rollers have to be contacted with the skin at same time to form a circuit.

Basic use method

1. Press the button for about 2 seconds to boot. The first indicator light is on.

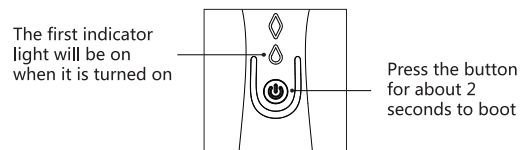

The indicator light indicates that the current nursing gear is from bottom to top 1 / 3


Click the key to select the gear position

1. Face plasticization and subcutaneous metal implantation should not be used; please use it under the guidance of a doctor during pregnancy.
2. The metal heads should be cleaned by paper after use to avoid the accumulation of grease and other residues. Do not wash it with water or contact with metal objects.
3. Do not put it in the bathroom and other wet places, avoid short circuit.
4. This instrument is suitable for the whole body, avoiding pustule type and inflammatory type acne.

Do not wash the fuselage and wipe it with a wet cloth



As shown in the following figure, move in the direction of the arrow in the order of $1 \leq 5$.

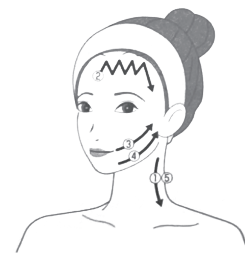

- 1.roll along the neck, from top to bottom
- 2.The middle of the forehead scrolls toward the temple direction curve
- 3.from the corner of the mouth to the back of the ear
- 4.roll along the contour line to the back of the ear
- 5.along the neck, scroll from top to bottom

As shown in the following figure, move in the direction of the arrow in the order of 1 / 3.



- 1.along the collarbone, moving from inside to outside
- 2.along the neck, moving from top to bottom
- 3.moving from the head-neck junction to the outside

- * EMS may stimulate dental nerves. Do not press the skin too hard and recommend the use of staggered areas. Especially when there are fillers in your teeth, you should pay more attention.
- * when used in the neck, note that the stimulation of EMS may be stronger.

Product name	Innovative Microcurrent 3D Roller
Product model	Amreetta
Product size	(host) L156*W62*H36mm
Product net weight	(host) 84g
Charging voltage	DC 5V / 0.3-0.4A
Current rated power	2W
Rated voltage / rated capacity	3.7V / 8 0 0 mAh
Charging time	about 3.5 hours
Strength gear	3 files
Production origin	China
Packing list	Roller Massager*1 bottom sitting*1
	manual*1 USB charging cable*1

hello@praimy.com