

正面

CULT

Vacuum Black Head Removal Device



USER MANUAL

www.primy.com



hello@primy.com

USER MANUAL

Blackhead Remover

Thank you for your purchase of our product. Please read the user manual carefully before using the blackhead remover for better enjoying the product.

Precautions

1. Blackhead remover is non-waterproof, therefore, please dry your hands before using to avoid device-body damages caused by entered liquid.
2. Please do not use other exfoliating products after using our device to avoid skin damage.
3. The usage frequency will be depended on your skin conditions. Please use no more than once per week for dry skins, and twice for normal and oily skins. Each time should be less than 5 mins.
4. Please do not press hard to make stronger suction, which might cause skin damage.
5. Please do not use the device for a specific zone constantly, which may cause skin purpura. Once the situation occurs, please stop using the device and contact us for consultation.
6. After using it, slight red skin is a normal situation. Usually, it takes about 30 mins to fade.
7. For the stubborn blackheads, please do not keep the probe sucking on one spot for too long. Such an operation may hurt the skin. Keep using this blackhead remover some more times can stop stubborn blackhead from getting back again.

ENGLISH

8. The device is not suitable for people who received cosmetic surgery or skin cosmetic surgery, and people with physical or mental disabilities. Children, pregnant women, and patients with heart disease should not use the device.

Quick Guide

Suction Head:
Using different suction heads to deal with different needs

Display screen:
Mode display
Low battery charging prompt


Button:
Short press to boot into the lowest power mode, and then press to upshift
Press upshift again

...
Press the button again at the highest mode to turn off the device.

Charging contacts:
Corresponding to the charging contacts of the charging base

ENGLISH

Charging indicator:
When charging, the back charge indicator flashing; When fully charged, the charge indicator lights up

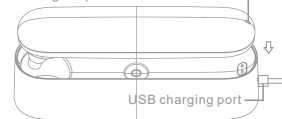

Suction Heads

Round suction head:
Features - Large suction power
Function - Removing Blackhead
Suitable for oily skin



Small round suction head:
Features - Small suction power
Function - Removing Blackhead
Suitable for thin and sensitive skin



Elliptical suction head:
Features - Medium suction power, special shape
Function - enhance tight skin
Applicable skin - all



Microcrystalline suction head:
Features - Small suction power
Function - Microdermabrasion, exfoliating
Applicable skin - all



ENGLISH