

Product Manual

FAVEOLA

Vacuum Deep Cleaner with HD camera Hot and Cold



PRAIMY



www.praimy.com

hello@praimy.com

Product Introduction

It is a multifunctional beauty instrument that can solve skin firming, hot compressing, cleansing pores, acne, blackheads and oils issues. We bring vacuum suction technology which can massage the skin and suck out the pores, blackheads, grease and dirt gently etc. At the same time, it can be lifted tightly to make the skin more tight and smooth, clear and shiny.

Vacuum Suction


According to the working principle of "vacuum negative pressure lifting technology", its unique vacuum suction system can lift the skin which can enhance the absorption of the skin so that the beauty liquid can penetrate deeply into the skin. The moisture makes the skin brighter and shiny. Meanwhile, it enhances the skin's defense and avoids dirt damage to keep skin tight and elastic.

Download the App

Download the app COOLEER from Appstore or scan the QR code to download the app and turn on the WI-FI in settings.



1

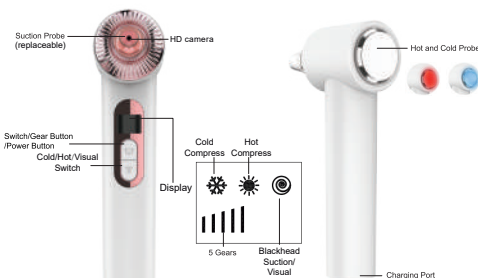

Usage Steps:



2

Beauty Probes

This product is equipped with six probes, each suction probe has different suction force, please kindly change the head from small to large, and gradually adapt to the suction force. The larger the probe is, the stronger the suction force will be, and the smaller the probe is, the weaker the suction force will be.

- Mini round hole suction head
Small suction, for small blackhead, sensitive skin, nasal groove acne, etc.
- Small round suction head
Small suction, suitable for soft and thin skin, less skin endurance.
- Medium and small round suction head
Medium suction, suitable for general oil secretion.
- Medium round suction head
Medium suction, suitable for acne with too much oil.
- Medium large round suction head
Large suction, suitable for the skin with strong endurance or more blackheads, need to move slowly on the skin after use.
- Large round suction head
It has strong suction. It is used to suck black head V face. When using it, it needs to move slowly on the skin. Do not suck more than 2S on one part. When using V face, it should be lifted from the bottom to the top and moved slowly from the chin to the top.

Safety Warning

1. Firming skin function has an effect of lifting, the skin which is hurt or has scars are not recommended to use.
2. Please do not press the tip to the skin, the product itself has a strong suction force.
3. After the nursing, due to the microcirculation of the skin, the skin will have a slight redness, which is a normal phenomenon and will disappear automatically after 30 minutes of inactivity.
4. Please do moisturizing after finishing the nursing.
5. Please do not use scrub or other similar products before or after treatment.

3

Function Usage Method

1. Download the app COOLEER from Appstore and turn on the WI-FI in settings. Press the "Switch/Gear Button" for 2 seconds. The product is turned on, the display is on, and the power is turned on. The suction display shows 1st gear. Press this button to turn the gear cycle mode from small to large.
2. Shortly press the "Switch/Gear Button" to turn on blackhead suction function. Short press plus one file, 5 files in total, in the fifth gear, another short press will close the blackhead function. Long-press the power button and switch to blackhead removal mode. Connect the WI-FI called cooleer, the password is 12345678, when the light inside the suction head is on, you can use this product.
3. Turn off the blackhead function. Shortly press the "Cold/Hot Switch" (M button) and the product enters into the cold compress mode. Gently apply the chilled probe to the area where the blackhead has been sucked to get the skin cold and help to shrink pores.

*Please turn off the product after using for 30 minutes or more.

1. Press "Cold/Hot Switch" again to switch to the hot mode, and gently apply the hot probe to the part to be sucked to help open the pores. Return back to the App see your skin condition when you're cleaning blackhead.
2. Adjust the blackhead probe to the appropriate position and slowly move the probe in the area where it needs to be treated. Do not suck the same part for a long time. (Before the tip leaves the skin, please open the suction head to the side to release the negative pressure air, then you can easily remove it.)

Product Parameters

Model: MZ-BR04	Gear: 5 gears
Net Weight: 202g	Working time: 120 minutes
Voltage: 5V 1A	Battery Capacity: 1000mAh
Power: 12W	Product Dimension: 190 * 45 * 45 mm

4

Precautions

- Make sure the filter sponge is properly placed in the filter cup before use, and the sponge should be kept dry.
- Do not wash the device by water.
- To avoid accidents, this product is not recommended to use for the following situations:
 - (1) Those who are being treated by a physician and being physically abnormal
 - (2) Those who have malignant tumor
 - (3) Pregnant women, maternity, children
 - (4) Infected patients
 - (5) Those who have skin sensation obstacle or abnormal skin
 - (6) Acute patients
- On the condition of the damage to the casing or mechanical or electrical failure of the equipment, please stop using this product.
- Please do not place the product in a place where the light is shining, and keep it in a dry and ventilated position.

Warranty Regulations

- The warranty of this product is 1 year, only for repair, maintenance for lifetime.
- During the warranty period, the failures caused by the normal use of the instructions (as judged by the company's official staff) can be repaired for free.
- During the warranty period, if one of the following conditions occurs, a repair fee will be charged:
 - (1) Can not provide this manual and valid purchase proof.
 - (2) Failure or damage caused by failure to use, disassemble, or send to the maintenance site.
 - (3) Failure or damage caused by improper storage, such as dropping, collision, and breakage after receiving the product.
 - (4) Failure or damage caused by other unavoidable external factors.
 - (5) Improper use of equipment causes water enters into the equipment or other damage caused by the solution.
 - (6) Damage caused by voltage using a power source instead of the designated one.
- If the product does not meet the warranty conditions and needs to be charged, the company will charge the consumer.
- This warranty card is not issued separately, please keep it in a safe place.

5