

# PURE

## Ultrasonic Skin Scrubber

[www.praimy.com](http://www.praimy.com)

[hello@praimy.com](mailto:hello@praimy.com)



USER MANUAL & WARRANTY INFORMATION



## Contents

01	.....	Introduction
02	.....	Important
03	.....	Important
04	.....	Product Mix
05	.....	Function Guide
06	.....	Function Guide
07	.....	Function Guide
08	.....	How to charging
09	.....	Cleaning the instrument
10	.....	Warranty terms and conditions, Environment
11	.....	Specifications, Disclaimer
12	.....	Appendix
13	.....	Appendix

## Introduction

Ultrasonic Skin Scrubber is a modern instrument that improves any skin care routine at home. It works like a Skin Scrubber to remove impurities from the pore of skin, with nutrition leading-in Function, and lifting skin function providing a radiant and flawless result without hurt. It provides unique functions and an exceptional experience for you.

## IMPORTANT

Read this user manual carefully before you use the instrument and keep it for future reference.

### Danger

- Keep the instrument dry. Never immerse the into water.
- Do not charge this instrument in wet surroundings (e.g. near a filled bath, a running shower or near a filled swimming pool).
- Never throw the instrument into fire.

### Warning

- Do not use the instrument, the attachment if it is damaged or broken to avoid injuries.
- Fully charge the instrument at least every 2 months to sustain battery lifetime.
- This instrument is not intended for use by persons younger than 12 years and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the instrument by a person responsible for their safety.
- Do not insert metal-containing material into the socket of the instrument plug to avoid short circuiting.
- Women who are pregnant or lactating are not recommended to use this instrument.


## IMPORTANT

### Caution

users to use products more securely and appropriately to prevent prevent potential accidents or hazards in use.

1. Please stop using this product if there are any side effects.
2. Patients with infectious diseases, heart failure, abnormal blood pressure and skin diseases should consult medical experts before using this product.
3. Please consult your dermatologist before using this product if you are using medical supplies or have been treated with abrasions.
4. Do not use this instrument if your skin is still recovering from a cosmetic procedure, such as a chemical or mechanical peeling, laser resurfacing, fractional laser treatment, or mesotherapy treatment or any other treatment with dermal fillers in the past two weeks.
5. Do not use in high temperature or humid environment.
6. Please keep this product away from water source and fire source; Please be careful to avoid electric shock during rainy days.
7. Charge, use and store the instrument at a temperature between 5°C/40°F and 35°C/95°F.
8. Do not use it for children or non-verbal communicators.
9. Do not drop, dismantle, repair or modify the product, as this will damage the pro-duct or cause a failure.
10. Before using this product for the first time, please charge it.
11. Do not use the instrument for any other purpose than described in this manual.

## Product Mix

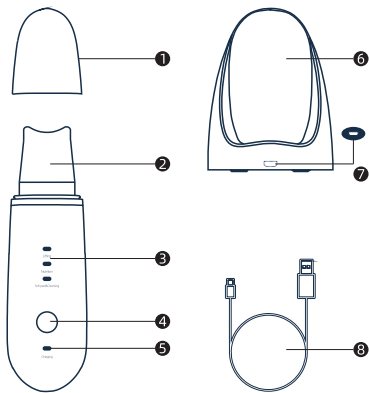

- ① Cover
- ② Spatula
- ③ Function Indicator Light
- ④ Power Switch(Switch Mode)
- ⑤ Charging Indicator
- ⑥ Charger
- ⑦ Charging Socket
- ⑧ USB Line

## Function Guide

1. It is composed of three modes: tile-tightening mode, nutrient import mode and deep cleansing mode.
2. Three models can be operated by a button, when the indicator light flashing for the pattern, click again, indicator will stop flashing and long on, represents in the pattern(powerful).  
**\*Auto Shut-off:** The instrument will automatically shut off after 10 minutes of continuous use. When using this instrument, if you wear gloves, it will result in no ion function.

### Lifting mode

1. In this mode, the positive ions produced by the instrument can reduce the pores, increase the elasticity of the skin, eliminate fine lines and wrinkles.
2. How to use

Step 1: press and hold the power switch for 3 seconds boot, and then choose the Lifting, when the indicator light flashing for the pattern, click again, indicator will stop flashing and long on, represents in the pattern (powerful).

Step 2: apply a nourishing cream or a facial mask, press the picture, and use a spatula to touch the skin and push the machine to massage and lift the skin from the bottom up.



## Function Guide

### Nutrition mode

1. In this mode, promote skin to absorb nutrients, such as vitamin C and min-erals, therefore, using this model, cooperate to use protect skin to taste the effect is more obvious, such as whitening lotion, contains vitamin C, vitamin A and vitamin E skin care products.

### 2. How to use

Step 1: press and hold the power switch for 3 seconds boot, and then choose the Nutrition, when the indicator light flashing for the pattern, click again, indicator will stop flashing and long on, represents in the pattern (powerful).

Step 2: apply a nourishing cream or a facial mask, press the picture, and use a spatula to touch the skin and massage.



It has the best effect on the skin circle, and it can massage 2-3s better.

## Function Guide

### Soft peel & Cleansing mode


1. The 25Khz ultrasonic wave generated in this mode is sprayed by water mole-cules, and the spray penetrates into the cuticle and is separated from the transparent layer.

### 2. How to use

Step 1: press and hold the power switch for 3 seconds boot, and then choose Soft peel & Cleansing, when the indicator light flashing for the pattern, click again, indicator will stop flashing and long on, represents in the pattern (powerful).

Step 2: misten the skin with gauze or cotton pad.


Step 3: press the picture, use a spatula to contact the skin, push the spatula and massage to clean the skin, a course of time is 10 minutes. After the course of treatment, the instrument automatically shut down.

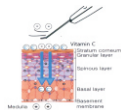

It has the best effect on the skin circle and push, and it can massage 2-3s better.

## How to charging

1. Connect the computer with the USB cable, turn on the computer, and connect the USB charging cable to the Charger, and place the ultrasonic skin scrubber above the charger , start charging.
2. Use the adapter to charge-connect the USB cable to the adapter, connect the USB charging cable to the Charger and place the ultrasonic skin scrubber above the charger , start charging.
3. The charging Indicator shines when charged, the charging Indicator is long bright when it is full of electricity.


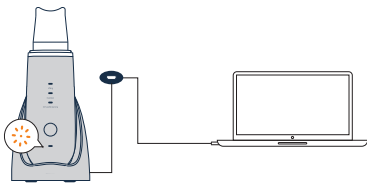


**CAUTION:** Before charging, make sure the charging socket and charger are free from water. Do NOT use your instrument while charging, and do NOT charge for more than 24 hours. Never use the instrument if it is damaged.

## Cleaning the instrument

Notes:

- Make sure the instrument is switched off before you start cleaning.
- Make sure the instrument is unplugged when you clean it. Only wipe the instrument with a dry cloth.
- Never use scouring pads, abrasive cleaning agents or aggressive liquids such as alcohol or acetone to clean the instrument.
- Clean the instrument and the Spatula thoroughly after each use to ensure optimum performance.
- Never rinse the whole instrument under water.
- Make sure all parts are dry before you reconnect to the adapter.

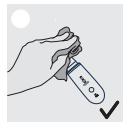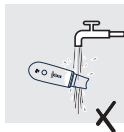

## Warranty terms and conditions

### 1-YEAR LIMITED WARRANTY

This instrument for a period of ONE (1) YEARS after the original date of purchase, against defects due to faulty workmanship or materials. If you discover a defect and notify us during the warranty period, we will, at its discretion, maintenance the instrument free of charge.

The warranty covers working parts that affect the instrument's function. It does NOT cover purely cosmetic deterioration caused by fair wear and tear or damage caused by accident, misuse or neglect. Any attempt to open or disassemble the instrument (or its accessories) will void the warranty.

For any warranty claim, reasonable evidence must be provided that the date of the claim is within the warranty period. For this reason, please keep your original purchase receipt for the duration of the warranty period.

## Environment

Do not throw away the instrument with the normal household waste at the end of its life, but hand it in and official collection point for recycling. By doing this, you help to preserve the environment.



This product must be sent to the official designated recycling center for recycling and recycling.

## Specifications

MATERIALS: ABS

SIZE: Instrument: 152\*46\*27 mm, Charger: 80\*61\*61 mm

WEIGHT: Instrument: 107g, Charger: 93g

BATTERY: Li-ion 600 mAh 3.7 V

CHARGING: 3 hours

USAGE: up to 100 min

STANDBY: 60 days

INTERFACE: 1-button

## Disclaimer

Users of this instrument do so at their own risk. Neither we nor its retailers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting, directly or indirectly, from the use of this instrument. Further, we reserves the right to revise this publication and to make changes from time to time in the contents thereof without obligation to notify any person of such revision or changes.

Please contact us for more comprehensive information, including user guide update, customer comments, frequently asked questions, etc.

