

Thanks for purchasing this product.


- Please read this manual carefully, and use this product safely and correctly.
- After reading, please keep the instruction manual for future reference.
- If this product was transferred to others, please be sure to attach this instruction manual at the same time.

Tips:
This product has heating function on the surface. For those who are not sensitive to heat, please pay attention that when using it for the first time, please test the temperature by hand first, adjust the temperature according to your own needs to avoid skin scald.

▼ Specification

Name	Face Lift Device with Gentle Brush
Model	VANESSA
Material	ABS,Silicone
Size	About 106*75*75mm
Weight	About 160g
Heating temperature	About 42°C ±6°C
Rated voltage	3.7V
Power	2.7W
Adaptor specification	5V/0.5A Type-C
Battery	Lithium battery
Charge time	About 1.5h

▼ Product structure

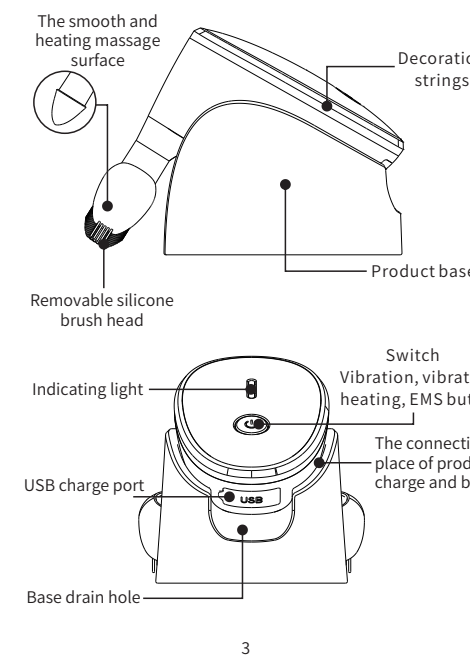

▼ Operation steps

The description of product button :

- 1: Long press the button to turn on the default mode: low vibration intensity (blue 1 light on).
- 2: Short press button for the first time, medium vibration intensity (blue 2 light on).
- 3: Press button second time, high vibration intensity (blue 3 light on).
- 4: Press button third time, heating function (red light on).
- 5: Press button fourth time, heating+vibration+EMS function (purple light one).
- 6: Press button fifth time, low EMS intensity gear (blue 1 light on).
- 7: Press button sixth time, Medium EMS intensity gear (blue 1,2 lights on).
- 8: Press button seventh time, high EMS intensity gear (blue 1,2,3 lights on).
- 9: Press button again, cycle above steps 1,2,3,4,5,6,7. Long press button to turn off.
- 10: Charge: Red light on when charge, blue light on when full charged.

▼ How to use

Before use
Please charge first: it can be connected to a computer or a power bank through a Type-C cable for charging.
Notice: Do not use the product when charging.

When cleaning the skin, use with a proper amount of cleanser, the effect will be better.
Note: When skin under dry and dehydrated condition, it is not recommended to use the product with power on, so as not to cause discomfort.

Charging operation method: Place the product on the base (the product USB port and the charging base are placed in the direction with Type-C), connect the data cable, plug in the power supply, and the indicator light is on.

▼ How to use

Product in use
(1) **Cleaning function**
Prepare the fully charged product, rub the cleanser or cleansing cream into the bubbles, apply it on the face and neck, turn on the product, set it to vibration mode, and stick the silicone heads on the area to be cleaned, drawing circles and making a circle way to clean, adjust the intensity of the gear according to your own situation.

(2) **Heating lead-in function**
Remove the silicone brush head, apply lotion or mask on the face, use the smooth massage surface to gently massage on the skin with a thermal gear for 3 minutes to open the pores and enhance absorption.

(3) **Vibration+heating+ems uplifting function**
Remove the silicone brush head, apply lotion or mask on the face, use the smooth massage surface to gently massage the

skin with heating, vibration, and EMS mode for 3 minutes to open the pores and enhance absorption.
(4) **EMS uplifting and tighten function**
Remove the silicone brush head, apply lotion or mask on the face, use the smooth massage surface to close to the skin, gently uplift from bottom to top, adjust the gear appropriately according to your needs.

Notice:
When using this product for the first time, please adjust it to the lowest gear. Even if you are very familiar with this product, before use, please adjust to the lowest gear and then increase the gear appropriately. Please do not use this product to achieve the purpose outside the scope of the description. It is recommended to use it for 3-5 minutes for each part.

The power of the product is too low, and the product will turn off after the charging indicator flashes for 5 seconds (to protect the battery from damage, if you do not use the product for more than 2 months, please fully charge it before placing it).

▼ How to use

After use
1. Long press to turn on/off the machine.
2. Use cloth to wipe the massage heads slightly to keep the product clean.
3. You can directly rinse the silicone brush head with clean water (please make sure it is turned off). After rinsing,

please dry it with a dry towel.
4. Although the product is waterproof, it cannot be immersed in water for a long time. The Type-C interface is not waterproof, please stay away from water.

▼ Common questions

1. The product cannot turn on. Please check if the product is fully charged. Whether the adapter has been removed from the power outlet.
It is recommended to fully charge the product before using it for the first time.
It is recommended to store it in full charged if not used for more than 2 months.
2. When using the function keys, no temperature is felt. Is the temperature gear button turned on?
3. Can I use the lead-in function without applying lotion and other skin care products?
In principle, it is not necessary to add lotion and mask, and it can also promote the blood circulation of the skin's subcutaneous tissue. But after applying the lotion or mask, the massage head can penetrate the nutrient essence into the skin's subcutaneous tissue better.
4. Is there a requirement for the using time and frequency of this product?
Use it at intervals of more than 24 hours a day, and each use time is 3-5 minutes. For the first 60 days, it is recommended to use it at least 5 times a week.

5. This product stops after using it for a while. Check if the product has electricity.
Note: Please contact the after-sales in time if there is an abnormal situation or the fault cannot be eliminated.

▼ Notice

- Although the product supports life waterproof, please close the silicone plug of type-C charging interface before use, and seriously damage the internal circuit of the product.
- Do not use clay bases, exfoliating creams, silicon-soluble granular cleansers, exfoliating creams or scrubs together with this product, so as not to damage the soft silicone brush head.
- Avoid direct sunlight, do not expose this product to high temperature environment or boiling water.
- For hygienic reasons, please do not share a product with multiple people.
- When the product is in working condition, please do not leave it alone.
- If there is any discomfort during use, please stop using it immediately.
- Keep away from fire and corrosive liquids.
- If the product is damaged, please do not continue to use it.
- It is strictly forbidden to disassemble this product by yourself.

▼ Warning

⚠ The following persons should not use this product

- People with allergies.
- Suffer from the following skin diseases or inflammations: Atopic dermatitis, Acne caused by sunburn, inflammation and other skin diseases.
- Patient with partial skin injury, inflammation or any other skin disease.
- Those with facial neurological disorders, those who cannot perceive heat.
- Pregnant women, or patients who have skin problems due to hormonal imbalances and are still undergoing treatment, please seek medical advice before use, and follow your doctor's advice.