

Thank you for purchasing this product! In order to use this product correctly, please read this manual before use; in order to read the manual at any time when you need it, please keep this manual in a safe place.

About this manual:
The product pictures and actual product features shown in this manual, as well as some rooms, can only be used as a reference.
The user will use the actual product as the standard, and the referenced data comes from the company's laboratory.

Safety Precautions
Be sure to read the safety precautions before use.
Please use this product correctly according to the contents and recommendations of the manual.
The product should be cleaned in time after use to avoid accumulation of grease and its residues.

Incorrect actions can lead to personal injury
Incorrect operation can result in personal safety or property damage

Danger

- Please use the charging cable supplied with this product or the charging cable that meets the product specifications, and connect the charger plug to the household socket correctly to avoid malfunction or fire caused by heat.
- Since this product has a built-in rechargeable battery, do not put it in a fire or use it in a high temperature environment, which may cause malfunction or fire due to heat.
- If the product is unusually hot or deformed, damaged, or smells, do not use it to avoid the risk of electric shock and fire.


Caution

- To prevent the occurrence of danger, do not use it with medical electronic equipment.
- People who have undergone heart surgery or carry medical electronic devices such as pacemakers or implanted defibrillators that are placed in the body, please consult your doctor before use.
- This product is not intended for use by persons with physical, sensory or mental problems, or those who lack experience and knowledge (including children), unless supervised by a guardian and instructed in the use of the product, and responsible for the safety of the user.
- If the machine is damaged, in order to avoid danger, it must be repaired by the manufacturer, its maintenance department or a specialist in a similar department.

- Do not use this product (such as alcohol, benzene, spray, etc.) near flammable or explosive materials, which may result in fire or explosion.
- Place packaging materials such as plastic bags where children can't reach them.
- Please turn off the power before cleaning and maintenance.
- Do not use this product for purposes other than those specified in this manual to avoid malfunction or fire due to heat.
- When using, if the product is abnormally hot, please stop using it immediately and shut it down.
- Do not use the machine if it is faulty.
- Do not immerse the appliance in water or rinse it directly with a water tap.
- If you feel any discomfort in your skin or body, stop using it immediately.
- Please do not store this product in a humid environment.
- When using, remove precious metal objects and watches.
- Please do not use it on skin parts with wounds.
- The following people do not use it, as this may cause accidents or skin or physical problems.
 - Metal allergies
 - Skin sensory impairment
 - skin infections, dermatitis, etc.
 - bleeding disorders (hemophilia, etc.)
 - People with heart disease, malignant tumors, epilepsy, diabetes, infectious diseases, acute conditions.
 - Long-term use of steroid hormone drugs or liver dysfunction caused by capillary vasodilatation.
 - In the early stages of unstable pregnancy, just after birth or during pregnancy, physical discomfort (with pregnancy induced hypertension syndrome, etc.)
- The following people please consult a doctor before using this product, failure to do so may result in an accident or skin or physical problems.
 - Skin sensitive
 - Being treated in a medical institution
 - Just undergoing surgery
 - Impaired in the body
 - In the past 3 months, surgery has been performed at the site of use, with long-lasting laser or deep chemical strapping.

If you are uncomfortable with the use of this product other than the above situations, please consult your doctor.

Part Identification

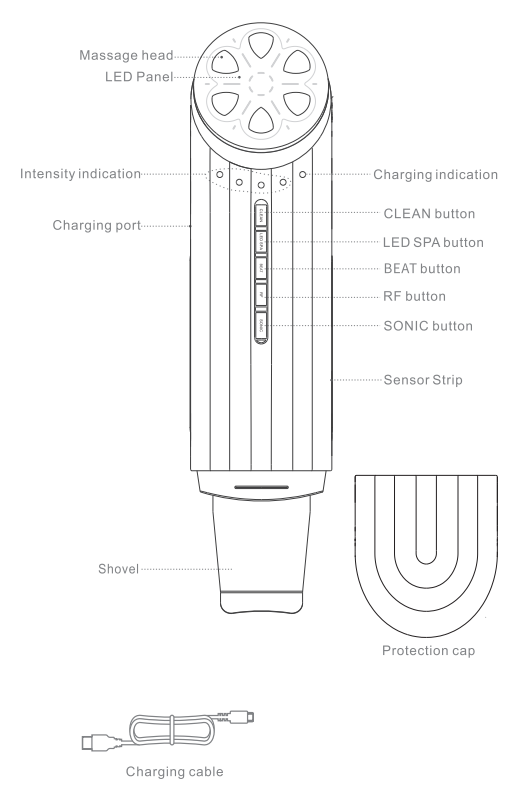

Before usage

Button function • Instructions

Strength/charging indicator

- From left to right, 1 light is on for first level, 2 lights for second level, 3 lights for third level, and 4 lights for fourth level. The higher the level, the stronger the output.
- The far right is the charging indicator (charging—blinking red light, charging completed—long bright red light).

Ion cleansing

Ionic current lead-out—deep skin cleansing.

LED skin care

Red LED—tighten the skin.
Can be used together with EMS and RF functions.

BEAT

Intermediate frequency intermittent electric pulse

Promotes skin nutrient absorption;
Can be used together with LED SPA function.

RF

RF firming

Make the skin feel warm and firm;
Can be used together with LED SPA function.

SONIC

Ultrasonic skin care

Clean skin—the default state when this function is turned on;
Nutritional skin care—After the function is turned on, short press again to switch to the flappping state.

“CLEAN”, “LED SPA”, “BEAT”, “RF”, “SONIC” each function can be turned on and off independently.

Intelligent • Timing

After the beauty instrument is turned on for 10 minutes, it will automatically turn off (there will be three beeps before turning off) and enter the standby state.
When operating each button of the beauty instrument, there will be a “di” prompt and a prompt that the LED light is on at the corresponding button.

Charging

When using it for the first time (for the first time or not used for a long time), please make sure to charge it first. And confirm that it has been fully charged (The rightmost indicator light is green) before use.


Charging time & times of use

Charging time	Times of use
About 2.5 hours	About 60 minutes

Please keep the charging environment at 25°C-35°C.
This is the operating time to be used immediately after full charging.
The service time will vary with the use mode, intensity and environment.
The charging time will be different due to the residual power of the charging battery and the use environment.

Charging method

- Insert the one socket of the charging line into the charging port of the device and the other socket into the USB port of the charger.
- Plug the charger into the household power socket.

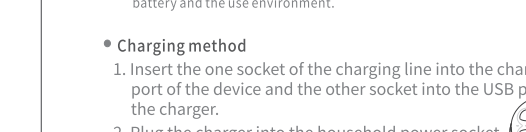

- Charging red light is shining
- Fully charged, red light is always on

- After charging is finished, pull out the charger from the socket, and then pull out the socket of charging line from the main unit.

This product only uses the provided charging line or the charging line and charger that meet the product specification requirements.
When using, plug the charger directly into the socket.
No extension line is allowed to be used.

Please be sure to confirm the plug shape and voltage (charger output is DC5V) before use, otherwise it may cause body heating, fire or failure.

Basic usage

[Preparation before Use]


- Please remove metal objects such as rings, necklaces, earrings, bracelets and so on.
- Please remove the makeup and clean the face before using it to ensure that the skin is in good condition.
- When performing eye care, please remove contact lenses first.



Skin care tips


- This product needs to be used together with beauty solution or beauty gel.
- Please do not use this product in larynx, eyeball, mucous membrane and so on.
- During use, do not stop in one place, or squeeze the skin hard.
- For the first time, it is recommended to start from the 1 gear and choose the suitable gear step by step according to the individual's skin feeling.
- Each function combination can be used together, please refer to the following nursing steps:
Ion cleansing → ultrasonic skin care → BEAT promotes skin nutrient absorption/red light can be superimposed LED skin care → RF radio frequency tightening (red light LED skin care can be superimposed).
- It is also possible to use a certain care mode separately according to personal needs, for example:
For those who often wear makeup, use ion cleansing and ultrasonic care after washing their face.
If skin care products have poor absorption effect, use RF to tighten after applying skin care products.
If you need to dilute the double chin, lift and thin your face, use BEAT to promote skin nutrient absorption.

[Ion cleansing]

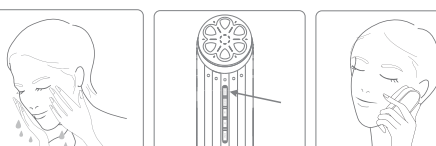


- Apply makeup remover or water to the skin.
- Press the “CLEAN” button to turn on ion cleansing.
- Slowly move the massage head from inside to outside and from bottom to top to clean the face.

[Ultrasonic skin care]

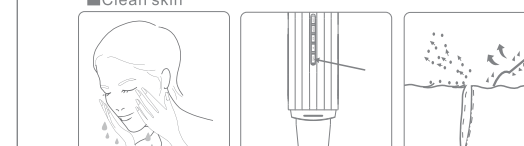

- Apply makeup remover or water to the skin.
- Press the “SONIC” button to turn on ultrasonic cleansing.
- Hold the machine with your hand, and slowly move the shovel blade (with the best position facing down) along the face to be cleaned.
- After the spatula comes in contact with the skin, water splashes (if the skin surface is not wet enough, you need to add water or cleansing lotion), the dirt and oil on the skin will be discharged out of the skin with the water.
- After cleaning the skin, press the “SONIC” button again to turn off the machine.

Nutritional skin care

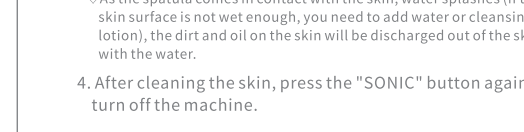

- Take the nutrient solution and apply it to the cleaned part.
- Press the “SONIC” button twice to turn on ultrasonic nutrition skin care.
- Reverse the shovel blade (the folding position of the shovel blade is downward), hold the machine with your hand and apply the shovel blade to the skin.
- After the treatment is completed, press the “SONIC” button again to shut down.

[Promotes skin nutrient absorption] BEAT 1-2 times a week

- Press the “BEAT” button to turn on the intermediate frequency intermittent electric pulse mode.
- After cleaning, apply the gel. If the skin is dry in the middle of the operation, continue to apply the gel or beauty liquid.

You can adjust the suitable gear by pressing the “BEAT” button according to your skin's feeling.
The LED skin care mode can be turned on at the same time, and red light skin care can be superimposed.

Hold the beauty instrument and slide it from under the chin to behind the ear, then massage to the shoulder. Then gently lift the neck upward.
Slide along the facial curve from the chin to the temples, then from the corner of the eye to the end of the eye to the temple.
Slide your eyebrow to the temple for a lift massage.




Once a day

[RF firming] RF

- Press the “RF” button to turn on the radio frequency skin tightening mode.
- After cleaning, apply the gel. If the skin is dry in the middle of the operation, continue to apply the gel or beauty liquid.

You can adjust the appropriate gear by pressing the “RF” button according to the feeling of your skin.
The LED skin care mode can be turned on at the same time, and red light skin care can be superimposed.

Place the beauty instrument in the in sequence for no more than 5 seconds, and then slowly slide to the next part.
Steps A to F are performed 2 to 3 times for each part.

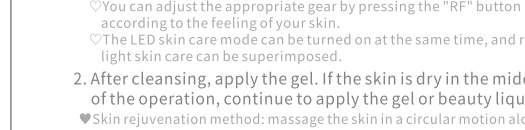

Cleaning and maintenance

- Please be sure to clean and maintain the product after each use.
- Cleaning or maintenance shall be performed under the shutdown state.
- Do not use steel wool, abrasive cleaner or diluent, gasoline, alcohol to wipe, otherwise it will cause scratch, deformation, rupture and other faults of components.

How to care for this product when the product is dirty

Wipe off stains such as finger prints with the cleaning cloth.

How to care for this product when the product is dirty

- Clean off any dirt. Rinse the product with cold or lukewarm water (35°C/95°F or below). If the product needs further cleaning, wipe with a wrung-out cloth that has been soaked in water-diluted detergent. After wiping, wash with water.
- Wipe off water droplets. Use a soft cloth to avoid scratching the surface of the product. Do not use an abrasive cloth such as metal polishing cloths, or an anything containing alcohol such as wet wipes.
- Air dry in a well-ventilated and shaded location.

Storage

- When storing the product, avoid a hot or humid location (such as inside a car or bathroom) and keep it away from direct sunlight.
- When not in use, keep in a safe place to avoid tripping on it and getting injured.
- Keep out of the reach of small children and store away from pets.
- Do not store near fire.
- Before using the product after a long period of disuse, check and make sure there are no problems with the product.
- If it is not used for a long time (more than 2 months), please take the product out every two months and store it after it is fully charged.

After-sales

- This appliance is intended for personal use only and is not suitable for commercial use.
- We offer a one-year warranty on the product, which is calculated from the date of purchase. During the warranty period, we will replace parts or the entire product due to the material or the production process.
- To apply for a warranty, the product must be shipped to an authorized repair station or factory along with a certificate (purchase invoice) that proves the date of purchase.
- The company does not guarantee the following conditions:
 - Damage caused by improper use, such as using the wrong power supply, inserting the product into an improper socket, man-made damage, etc.
 - The product caused by normal use is aging and wear, but it has little effect on the function of the product.
 - Product problems caused by other uses, or improper handling or non-professional repairs or the use of accessories other than the company.
 - Consumables or bonus accessories.
 - Damage caused by force majeure factors (such as lightning, earthquakes, fires, floods, etc.).

Waste Treatment

The marking on product or packaging material or accessory indicates that the material is made from a recyclable material. Separate waste into different recycling levels to facilitate the recycling of useful materials.

There is a lithium battery inside the appliance, therefore:

- The battery must be removed from the appliance before it:
 - When removing the battery, the appliance must be turned off;
 - The battery should be disposed of safely;
 - Users cannot replace the battery by themselves.
- Steps to remove the battery:
 - Pry open the front and back cover of the product from the assembly seam;
 - Unplug the battery leads and remove the battery.

Product Name: Multifunctional Premium Skin Scrubber
Rated Power: 5.0W
Net Weight: 127g
Charging Power

Specification

Product Name: Multifunctional Premium Skin Scrubber
Model: FANCY
Rated Power: 5.0 W
Net Weight: 127 g
Charging Power Supply: 5V—1A
Packing List: 1 main unit, 1 charging cable, 1 operation manual.


Due to the improvement of the products, some specifications may be changed without prior notice.

Fault analysis and elimination

In case of any malfunction during the operation, please carry out following inspections. If failures cannot be eliminated or the product cannot operate properly, please contact after-sales service center. Do not disassemble and attempt to repair it by yourself.

Situation	Possible Reason	Measure
Charging light is not lit	Battery is dormant The charging wire or charger has been damaged.	The battery needs 6-30 minutes back to working. It is a normal phenomenon and needs to be fully charged before use. The battery needs 6-30 minutes back to working. It is a normal phenomenon and needs to be fully charged before use.
Can't turn on the power after it off after a "di" beep	Low battery	Connect the beauty instrument to the charger and fully charge it before using it.

FANCY Multifunctional Premium Skin Scrubber

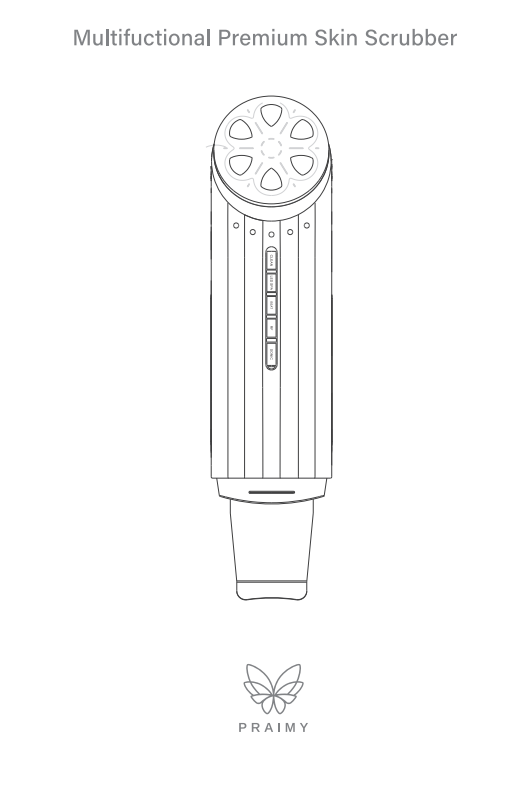

he@praimy.com