

BLAST

5 in 1 Face Lift Combination

USER MANUAL

Edition: 20210722

1

5 in 1 Rechargeable Ultra Skin Rejuvenator

Let Everyone Enjoy suitable home beauty device

! Notes

- Do not exceed the scope of use specified in this user manual;
- Pregnant women should follow doctor's advice;
- Do not use near area with metal prosthesis, such as eye balls, inner ear, bone areas and joints;
- Prohibit use by patients with pacemaker, osteoporosis patients, without sense of temperature or pain, and patients with skin disease;
- Prohibit use in vein or arterial pathology areas or Patients with Thrombophlebitis or Blood Coagulation Precedents;
- Prohibit use for patients with transmute changes, alcohol abuse and liver infections;
- Avoid physical impact;
- Only use according to the user manual; Prohibit disassemble, transform, maintenance or retrofit the product;
- Please store the device in a clean and stable place; Do not store in a place overheated or over cold; Keep away from fire, such as stove;
- The device is not waterproof, keep it away from water source;
- Unplug the plug after charge;
- Keep the device away from children;
- Warranty 1 year upon delivery, under situation of normal use.

2

5 in 1 Rechargeable Ultra Skin Rejuvenator

Let Everyone Enjoy suitable home beauty device

Instruction

Thank you for purchasing our product. To fully benefit from the support that our offers, please read the user manual carefully before you use the appliance and keep it for future reference.

Specification

N.W.:165g
 Size: 190(L) x 60(W) x 50(H)mm
 Working: environment: 5-40 °C
 Power supply: AC100-240V 50/60HZ;
 DC5V $\overline{\text{---}}$
 Ultrasonic frequency: 3MHZ
 LED RED: 625±5nm
 LED BLUE: 465±5nm
 LED GREEN: 525±5nm
 Charging time: Approximately 2.5H

Full set including

1 x Main device
 1 x USB Charging cable
 1 x Pouch
 1 x User manual
 1 x Gift box
 1 x cotton holder


2

3

5 in 1 Rechargeable Ultra Skin Rejuvenator

Let Everyone Enjoy suitable home beauty device

Button and accessory instruction

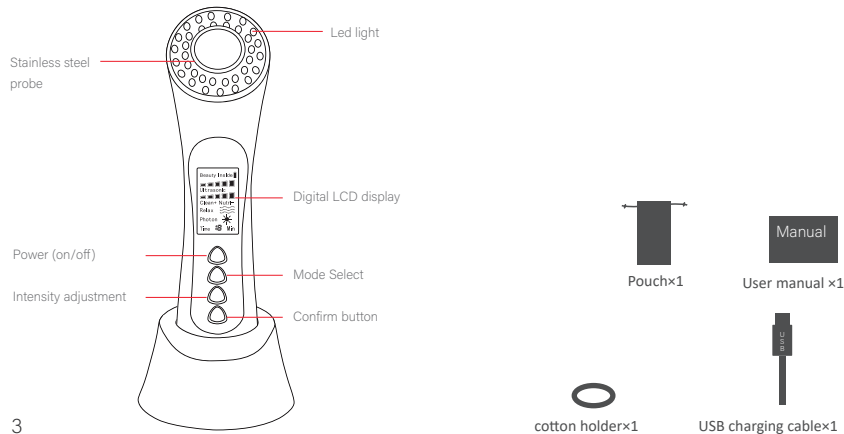

3

Operation instruction



1. Power On/Off the device.

2. Mode selection: Press this button to select Mode. The related mode will flash when choosed. Then press Intensity Adjustment button to choose intensity.



3. Intensity adjustment: When press Mode Selection button and related function is flashing, press this button to choose intensity. 1-5 levels adjustable. For vibration two modes for adjust, and time adjust with 1-15 mins.



4. Confirm button: After finish select the functions you need, press Confirm Button and device starts working. You can press this button to pause and adjust levels or choose other functions again.



Function instruction

1. 3MHZ Ultrasonic Function: Ultrasound is a high-level sound wave technology that penetrates deeper than the dermal layer. In the skin it increases the temperature of the surrounding tissue, thus increasing oxygenation to the cells, stimulating collagen protein and aiding in the removal of waste products. This is an anti bacterial cleansing action. It benefits the skin by stimulating blood circulation, diminishing fine lines and wrinkles, clarifies acne and blemishes and can also eradicate irregular pigmentation and redness.(Portable Galvanic ultrasonic beauty & health instrument Erasing wrinkles).

2. Galvanic function: Positive electrode can generate acid action with the functions of tranquilizing nerves, dropping blood supply in skin, strengthen the fiber tissue of skin, shrinking pores and relieving reddish skin. It is applicable for allergic skin and acne.

3. Galvanic function: Negative electrode can generate alkali action, stimulate nerves, promote blood circulation, soften fiber tissue and strengthen firmness, and is applicable for dry and ageing skin with slow metabolism.

4. Vibration Function: Microwave frequency vibration 12,000 rpm deep dermis, can soften the residue of cosmetics in pores, grease, blackheads, help to facilitate the discharge of oil secretion, deep export acne and blackheads, pimples and acne, and the formation of cosmetics residues, to prevent the formation of acne blackheads, spots.

6

5 in 1 Rechargeable Ultra Skin Rejuvenator

Let Everyone Enjoy suitable home beauty device

Different functions for different lights

5. Bio Red Light: wavelength $625 \pm 3\text{nm}$: Whitening and lighten spots, tendering skin and dispel crinkles, repair damaged skin, heal slender wrinkles, tighten pore, and proliferate collagen protein.

6. Bio Blue Light: wavelength $465 \pm 3\text{nm}$: diminish inflammation, kill bacterium and propionibacterium acnes, remove acne and scars etc.

7. Bio Green Lights: wavelength $525+/-3\text{nm}$: has the effects of neutralization, balancing and calming. It can ease strain, dredge lymph and useful in edema elimination.

7

5 in 1 Rechargeable Ultra Skin Rejuvenator

Let Everyone Enjoy suitable home beauty device

Skin care steps

Daily skin cleansing before using the device, apply your lotion/cream on your skin then use the device as following below picture, the probe fully should contact with your skin. Please avoid direct light to your eyes. Clean the probe with wet wipe and storage the device in a dry & cool place after use.

The recommended treatment time:

At most three times per week, no more than 15 minutes each time, to achieve the best results.

Note: Move the device gently over your skin during treatment.


Notice:

During the treatment, DO NOT look directly at the light-waves TO AVOID discomfort or dizziness.



How to Charge

The device will be charged and discharged for a few times before delivery.
Battery left about 50% after transportation.
Please full charge the device before use.
Connect to power supply for charging.

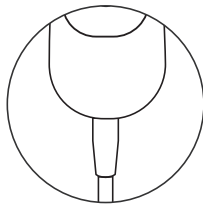

Maintenance and storage

- Cleanse & Dry the device after each use;
- Storage the device in a dry & cool place after cleanse;
- It is not recommend to share the device with others;
- Keep the device and heads dry and clean after each use;
- Keep this user manual for reference;

Malfunction & maintenance

If the device does not work properly, please check following :

- Make sure the device fully charged.
- Make sure press the "Power" button and Working indication light lights up;
- If no reaction after power on the device, maybe low battery; Please full charge it first and then try again;



www.praimy.com

hello@praimy.com

