

# Instruction Manual




## MOON EMS Neck Lifting Massager



www.primy.com

hello@primy.com

### Product Structure

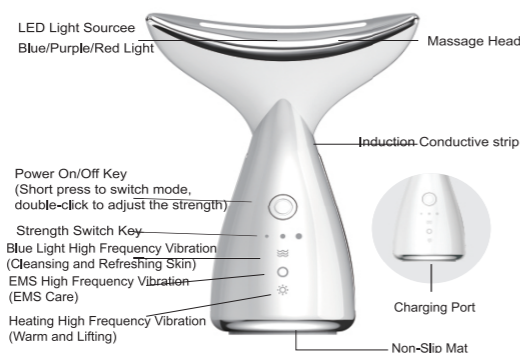

1. Long press to turn on/off, short press to switch mode, double-click to adjust the strength
2. Each mode is scheduled for five minutes, if it scheduled for 4 minutes, the corresponding mode indicator will flash.
3. When the battery is below 3.3V, the three mode indicators automatically flash 10 times and power off automatically.
4. The device is not working when charging, and the intensity indicator lights are flashing gradually, the three lights are often on after full charged.

Name: EMS Neck Lifting Massager	3.7V Rated Voltage: 3.7V
DS-8822 Model : MOON	Material: ABS + Plating
Unit Size: 85*106*45MM	1.5W Rated Input: 1.5W
3.5H Charging Time: 3.5hrs	Battery Capacity: 500mAh
Main Functions: Lifting and Tightening, Invigorating Skin	

### Mode Declaration

Cleansing and Refreshing Skin Blue Light + High Frequency Vibration	EMS Care Purple Light + High Frequency Vibration	Warm and Lifting Constant Temperature 45°C + High Frequency Vibration
--	---	---

- Blue light** Blue light can promote the synthesis of collagen, activate the skin, tighten the loose skin, collect fine pores, and have a great improvement to the oily skin and serious sensitive skin.
- Red light** Red light can enhance cell vitality, accelerate blood circulation, promote the growth of fibroblast and collagen, make the effect on smooth fine lines, wrinkles, tighten the skin, reduce light pigment, desalinate the scars and so on.
- Purple light** Purple light can improve cell oxygen, use alternating function to promote body surface microcirculation, drag the lymph and drain the edema, and have the effect of calming the skin, balance and stability.
- EMS** EMS micro-current releases intermittent pulse current, stimulates ATR production, improves facial muscle fiber elasticity, activates collagen activity, and tightens the skin.

### Product Features

The streamline massage head is specially designed for the neck, and the streamline design is in line with human industrial theory and transfers energy to every inch of the neck skin.

Massage head sonic vibration with three modes, accompany with blue, red, purple three colors light, EMS massage, jointly promote blood circulation, stimulate the acupoints of the facial neck, dredge meridians, improve circulation and promote metabolism, as to achieve desalination of neck lines, reduce edema, relax shoulder and neck, acupoint therapy, quickly eliminate fatigue and restore physical strength.

The three modes are: tendering skin and anti-wrinkle (blue light), warm and lifting (red light), EMS tightening (purple light), long press the power on/off key to start the device, convert the mode by each press, until a cycle is completed, long press the power on/off key to shut off.

There are three intensities with each mode, and each mode comes with the vibration.

USB cable charging, safe and environmental protection.  
Mini portable in shape, for use anytime and anywhere.  
The new introduction of the bionic design concept, fish tail shape design, delicate and novel.

#### Notes

Full charged before using the device.

If the product is not used for a long time (about 3 months), please take it out and full charged to store it.

### Neck Care Steps

Apply the skin care products evenly to the neck

Choose the mode according to personal needs, the massage head is pulled and lift from the bottom to top in the neck.

The back side neck also needs to massage, massage with press, pull, caress, push and other ways in the acupoints.

### Face Care Steps

Apply the skin care products evenly to the face

Choose the mode according to personal needs

Acupoints press to relieve muscle tension

Fit the facial skin, pull and lift from the bottom to top, make the face in shape.

### Attentions

Short press the power on/off key 1 second to shut off the device

Wipe the device with cleaning cloth and do not use corrosive chemical cleaning agent.

Place it in a dry and ventilation area.

### Attentions

Before using the product, please read the instructions carefully and keep it properly for future reference.

The product is charged with USB line, which takes about 3.5 hours to full charged when the power is insufficient, and can be used for 2-12 hours continuously.

Flash red light when low power, and it display red light when charging, green light after full charged.

The product has memory mode function and automatically enters the mode of the last shutdown once restart.

This product is not allowed to use when charging.

The device has a heating function, about 45 degrees, and the skin has a little difference in the temperature value, it may feel hot, after the skin applies the essential oil/massage cream, the device will automatically adjust to the most comfortable temperature.

If you feel unwell during use or have abnormal skin, please stop to use it and consult your doctor immediately.

The device has no waterproof function, do not place it in water.

Please place the device in a dry and ventilation area to avoid direct sunlight, free from corrosive chemicals.

When taking beauty care, please take the right amount of skin care products. If excessive use of skin care products, the skin cannot be fully absorbed and lead to waste, and there is the possibility of residual fluid infiltrating into the interior of the product, make the damage to the service life of the product. Please wipe it with a clean cloth after use.

Do not permit the children touch this product to avoid danger.

### Common use Q & A

Questions	Solutions
The device suddenly did not vibration, did not work, what happened?	Please check if the device power is low. Is there a correct boot on? After charging, restart the product. Please contact the brand dealer for consultation, or call the service hotline for after-sales consultation.
May I use the device every day?	Yes, you can. The vibration strength and temperature adjustment of this device is designed by a professional engineering design team, which can be used many times a day.
Allergic skin, can it be used?	Especially choose environmental protection, conform to international standards of materials, safety, health, no irritation and allergy to the skin, please rest assured to daily use.
At what time is it used every day?	Available for any time period, recommended to use it 5-10 minutes every time.

### Certificate

Date of Production \_\_\_\_\_

Examination Clerk \_\_\_\_\_

Operative Norm. GB4706.1-2005