

TWIST
MICROCURENT LED THERAPY NECK MASSAGER



Directory

Directory	01
Introduction	02
Product structure	03
Method of use	04
Method of use	05
Functional description	06
Introduction of Color Light Therapy and Vibration Therapy	07
Usage frequency and cycle / cleaning and maintenance use and storage considerations	08
Taboo population	09
Frequently asked questions	10
Product specification	11
Environmental conditions requirements	12
Repair regulations	13
guarantee card	14


Introduction

V face beauty instrument is a portable beauty salon equivalent care skin care instrument, using six modes of color light therapy deep into the true layer of the skin, stimulate cell growth and blood circulation, replenish moisture to the skin, delay skin transformation Combined with vibration mode, tighten the outline, improve the three-dimensional feeling, shape the V face, at the same time improve the muscle strength of fine muscles, enhance blood circulation, help to repair and maintain elastic skin, so that the outer layer of the skin nutrients infiltrate into the skin dermis. Thin face at the same time let facial skin enjoy different SPA massage.

Before using the V face beauty instrument for the first time, please carefully read and understand the instructions, instructions, and safety information of this nursing instrument, and use it in accordance with its guidelines.


If any query, please contact after sale.

Product structure

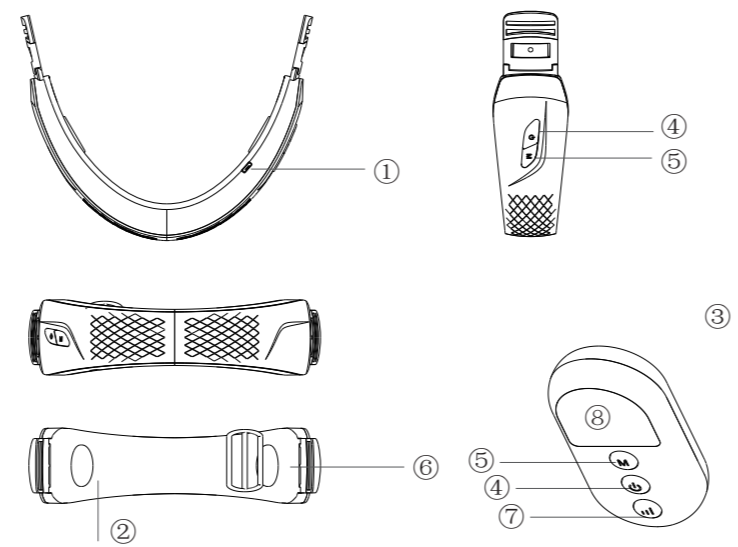


- ① Charging Port
- ② Intense Pulsed Light (IPL) Emitting Area
- ③ Remote Control
- ④ Power On/Off
- ⑤ Mode Button
- ⑥ Charging Indicator
- ⑦ Strength Selection
- ⑧ Remote Control Display


Method

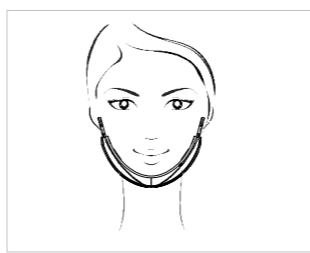

1 The V face beauty instrument can be controlled by a wireless remote control device or use alone, you can go on follow-up operation after fitting the instrument to your face, and there is a telescopic device at the bottom of the instrument (the place to be in touched with the chin), which can be adjusted to fit the shape of your face perfectly.

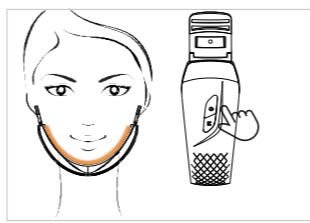


2 Long press power button on instrument or the remote control for 1s to turn on, the instrument can be used independently. Press the instrument or remote control, the buzzer beeps once for turning on, and turning off if beeps twice.

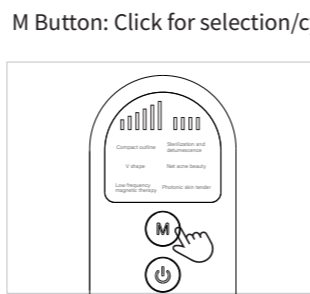


M Button: Click for selection/cycle
6 Modes Available In Succession: Contour Tightening, V-shaped Face Shaping, Low Frequency Reduction, Acne Cleansing & Repair, and Photorejuvenation. You can press button "M" to adjust to the next mode.


Method

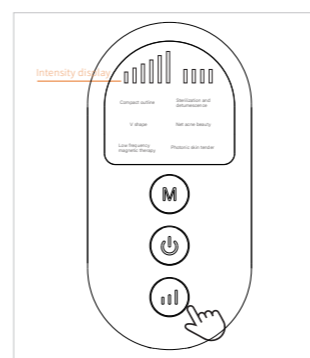

4 The intensity can be adjusted by intensity Button (). The default intensity of the device is "NONE". Each press to the intensity button, an intensity will be increased. One more press can return to "NONE" after reaching the highest intensity (long press to cancel all the current intensity), then you will start the second cycle press.

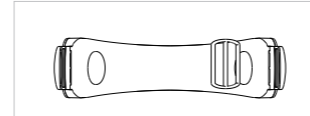

5 Take off the instrument and store it in a cool and dry place.

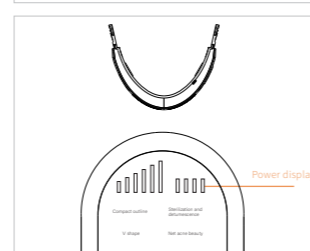

6 The instrument would stop working when charging. Redray on when charging, and Bluearay on when full charged. The battery symbol of the remote control: All 4 lights are on when full power, the first light would flash when the battery is in low power, and it would shut down when 5 times flashes occurred. If the device is exhausted, please plug the USB cable into the device and connect it to the power source for charging.

Function Description

- Contour Tightening: EMS+Bluearay
- V-shaped Face Shaping: EMS
- Low Frequency Magnetic Therapy: EMS+Redray
- Sterilization & Swelling Reduction: Redray+EMS+Heating
- Acne Cleansing & Repair: Redray
- Photorejuvenation: EMS+Redray+Vibration


Daily Care—Intermittent Alternating Redray & Bluearay and Vibration Mode: In this mode, firstly, the device will irradiate 3 times with Redray from left to right and right to left, and then switch to Bluearay irradiation for 3 times, so as to cyclically perform Redray and Bluearay irradiation vibration therapy, which can be used for daily skin cares.

Emergency Skin Repair—Continuous Alternating Red & Blue Light and Vibration Mode: In this mode, Redray and Bluearay are continuously switched from the middle to the outsides to irradiate and vibrate for therapy. It can be used for emergency repairing and cares, and it needs Bluearay and Redray quickly sync for skin repairing.

Smart Skin Care—Bluearay Flashing Irradiation Mode: In this mode, the Bluearay flashes slowly and illuminates the skin without vibration. It can be used for over-fatigued skins, reduce skin burden and soothe the skin.

Introduction to IPL and Vibration Therapy

Four major therapies: Redray Irradiation, Bluearay Irradiation, EMS Micro Current, Vibration. Use alone or in combination to enhance skin elasticity and shape a V face.

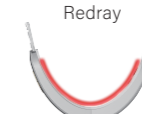

Redray Irradiation: Improve active of skin cells, promote metabolism, stimulate secrete a large amount of collagen and fill the fibrous tissue itself. Help promote blood circulation, enhance skin elasticity, anti-aging, anti-oxidation and skin repair.



Bluearay Irradiation: It can quickly inhibit common skin inflammation, calm the skin and regulate the secretion of keratinocytes, clean and sterilize the skin, strengthen skin absorption and repair, balance skin oils, whiten spots, rejuvenate wrinkles, and reduce pores.



EMS Micro Current: Release intermittent pulse current to stimulate the production of ATR, stimulate human subcutaneous tissue through EMS micro current, improve facial muscle fiber elasticity, activate collagen activity and tighten skin.



Vibration: Soothes facial muscles, improves the strength of fine muscles, increases facial elasticity, fully lifts and tightens facial muscles, help achieve shaping a V face.

Frequency and Period

Each mode is set for 15 minutes, and it will automatically stop working after 15 minutes, the previous use time would be adding up. If you need to continue, restart and selection mode.

Recommended to use the device 2-3 times a week regularly for better effect.

Cleaning & Maintenance

To better preserve the product, we recommended to clean the equipment regularly, especially areas touched with skin directly.

1. Turn off the power completely after using, and gently wipe the instrument with a clean dry cloth.
2. Please do not wipe and clean by flammable liquids such as alcohol (including perfume, disinfectant or other substances containing alcohol), using acetone or other chemicals.
3. When no usage over 7 days, separating the regulator and the charging cable, and be sure to keep it power off.

Cleaning & Maintenance

4. When low or no battery power, be sure to turn off the power before connect the charging cable for charging.

Please do not immerse the device or charging accessories in water.

Tips For Use And Storage

1. Forbidden to use besides the original purpose.
2. Excessive use may cause a burden to skin. If adverse reactions occur, please stop using.
3. Forbidden to use and store in high-temperature places such as kitchens or cars, dusty places, and places with fireworks, which will shorten the life of the product.
4. Charge the product in time when the power is insufficient to avoid affecting the product life.
5. Avoid multiple drops and strong collisions with sharp objects.
6. It is prohibited to disassemble or repair the product by yourself.
7. When storing the product, be careful to choose a place that children cannot touch.

Taboo People

- a. Excessively exposed skin, advised to wait for skin repair before using, so as not to cause secondary damage to the skin.
- b. There is obvious trauma or wound on the face.
- c. Patients with severe skin inflammation or sensitive skin are advised to have a 1 minute test first, and continue to use without adverse reactions.
- d. Facial fractures or major post-operations need to be used under the guidance of a doctor.
- e. Use with caution in bleeding areas.
- f. Patients with eardrum inflammation or hearing aids need to be used under the guidance of a doctor.
- g. Patients with metallic diseases.
- h. Patients with severe heart disease.
- i. Blood pressure disorder.
- j. Children under 12 years old.
- k. Pregnant women need to be used under the guidance of a doctor.
- l. Face cosmetic surgery needs to be used under the guidance of a doctor.

Q&A

- 1) Q: What should I do if the device does not work at all?
A: The battery is insufficient, please restart the product after charging is completed.
- 2) Q: What is the cause of the automatic power off during use?
A: Low battery. Please turn off the product and use it after charging.
- 3) Q: How long for each use? and how long to be effective?
A: According to the individual's facial condition, 2-3 times a week is advised and regular use would be better.
- 4) Q: Does the instrument have a memory function?
A: Each mode has a memory function. When the user adjusts to other gears, the previous gear will form a memory, and it will default to the last used frequency.

Product Specifications

Product Name	Microcurrent LED Therapy Neck Massager
Product Model	TWIST
Charging Voltage	DC3.7-5.0V

• Remote Control

Battery Size	35*70*15mm
--------------	------------

• V-face Lift Belt

Product Size	177.8*106.3*49mm
Charging Time	3 hours
Rated Power	3.5W

Environmental conditions requirements

Working Conditions

Temperature	+10 °C - 35 °C
Relative Humidity	3% - 75%rH
Air Pressure	700 - 1060hpA

Storage Conditions and Transport

Temperature	-40 °C - +70 °C
Relative humidity	10% - 90%rH
Air pressure	500 - 1060hpA

Repair provisions

This book promises to be free of charge for repairs within one year from the date of purchase. In compliance with the use of instructions and other matters needing attention in the normal use of the state of failure, according to the contents of the warranty will be given free repair. Please bring this product and warranty to your point of sale at the time of purchase or the local customer consultation center to apply for repair.

1. Use this product in compliance with the instructions, such as failure during warranty period, the company will give free repair (except accessories)
2. even during the warranty period but in principle the following cases are paid for repairs.
 - ① Failure and damage caused by misuse other than instruction manual.
 - ② Failure and damage resulting from improper repair and modification.
 - ③ Failure and damage caused by place of storage, transportation, drop, etc.
 - ④ Failure and damage caused by use outside the home (e.g. business).
 - ⑤ No instructions can be provided.
- ⑥ The warranty card does not provide the date of purchase, the name of the customer, the name of the shop, or the handwriting is altered (please provide the receipt of the invoice for the purchase of correspondence)
- ⑦ Failure and damage on board a vehicle.
- ⑧ Fire, earthquake, flood, thunder, other natural disasters collapse, abnormal voltage, use of power sources other than specified (voltage, weekly wave number), etc. Failure and damage caused by faulty repairs and modifications carried out outside our company.

This specification will not be reissued, please keep it properly.

Product warranty card

The user keeps it and shows it to the maintenance site during maintenance

Product Name:	
Type Number:	
Production Land:	
Duration of warranty:	
Invoice No.:	
Date of purchase:	

Maintenance records

Delivery date	Date of submission
Troubleshooting, maintenance, handling:	
Maintenance unit:	Signature of maintenance personnel:

www.praimy.com

hello@praimy.com