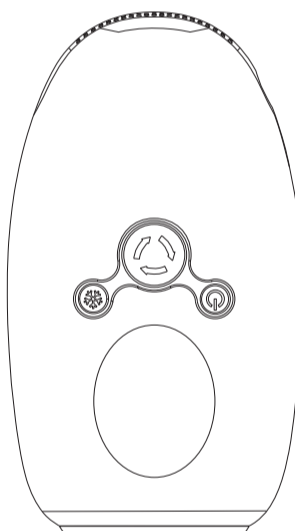



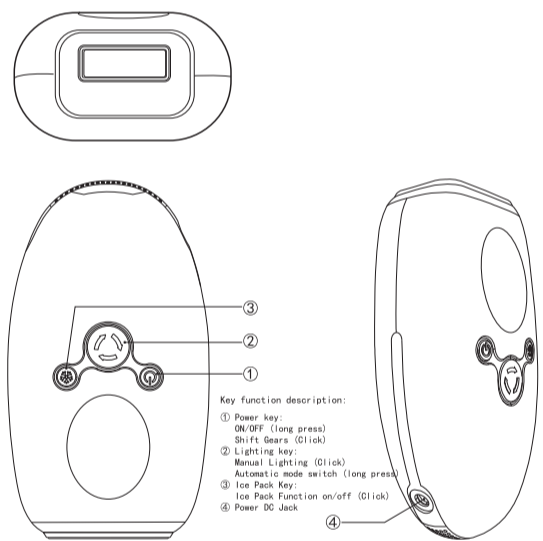
ZENANA
Ice Cooling IPL Hair Removal



PRODUCT INSTRUCTIONS

www.primy.com  hello@primy.com

1. Device Introduction



Device Name	Ice Cooling IPL Hair Removal	Device Size	136.5 x 77 x 47.5mm
Power Adapter	100-240V	Illumination Area	3cm ²
Rated Input	DC12V/3A	Lifespan	990000 times
Max Output	36W	Net Weight	Around 225g

2. Working Principle

- Under the illumination of impulse light with specific wave band and pulse width, the melanin in the skin absorbs energy and transfers the transient heat to the follicle. When the hair follicles reach certain temperature, they will enter the dormancy period, which makes the hair stop growing by stages and fall off. The ideal effect can be achieved by using the remover for 8 weeks or more, decided according to the physique of different people. The built-in filter will filter out ultraviolet light and other useless light, which makes the device ideal for skin care.
- Ice-cool function can efficiently relieve the pain (any feeling may vary), thus improving the experience of our product.
- To achieve the best effect, please keep using this device for around 8 weeks; At that point, hair would be thinner and start to fall off. After that, remember to reuse it once within 15-30 days to maintain the effect depending on how fast the hair grows back.

(Note: even for the same power level, burning sensation may vary from person to person)


3. Suitable skin Color

<input checked="" type="checkbox"/> White or alike	<input checked="" type="checkbox"/> Beige
<input checked="" type="checkbox"/> Light brown	<input checked="" type="checkbox"/> Medium brown
<input checked="" type="checkbox"/> Dark brown	<input checked="" type="checkbox"/> Dark or alike

4. Must know

- Do not charge the device since it's plug in. Please unplug the power cord after usage and well store them in the shade.
- Please shave the hair off first before usage and please go with shaving cream. Otherwise, the skin would be damaged under dry condition.
- Due to different kinds of skin, for first time user, please run a trial test on small area then use it on a large scale if without adverse reactions in 12 hours.

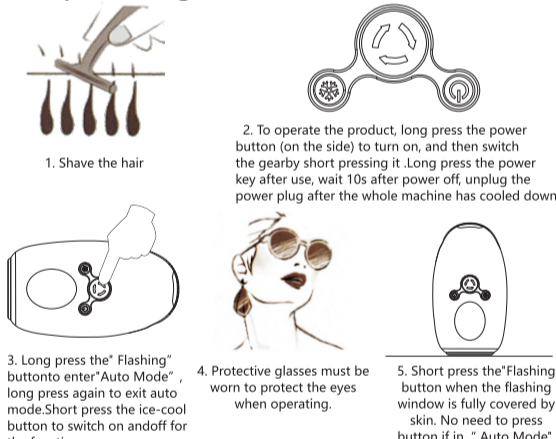5. Operating Guide






- Shave the hair
- To operate the product, long press the power button (on the side) to turn on, and then adjust the gear (short pressing it). Long press the power key after use, wait 10s after power off, unplug the power plug after the whole machine has cooled down.
- Long press the "Flashing" button to enter "Auto Mode", long press again to exit auto mode. Short press the ice-cool button to switch on and off for the function.
- Protective glasses must be worn to protect the eyes when operating.
- Short press the "Flashing" button when the flashing window is fully covered by skin. No need to press button if in "Auto Mode".

6. Not Applicable

- Female under pregnancy, period and lactation, light-sensitive person, juvenile and people with cicatrix are not allowed to use.
- People with heart disease, diabetes, hypertension, malignant tumor, folliculitis, infectious skin disease, burned skin and anticoagulation disorder are prohibited.
- If other methods of hair removal are applied lately, please do not use this device until hair grows back.
- If products with whitening element or Vitamin A are applied, please wait one week from the day of the application.
- This device can be used with essential oil products.
- Patients taking aspirin, immunosuppressive drugs, oversensitive drugs (such as tetracycline, sulfonamides, etc.) are not recommended to use this device for the time being.

Hair removal cycle

Week 1-4	Twice a week	Hair grows slower	
Week 5-8	Once a week	Hair becomes thinner and softer and some of the hair start to fall off	
Week 9 or above	Once within 15 to 30 days	The hairs start to shrink and hair stops growing	

7. Attention

- People with insolation or sunburn within last two weeks must wait until the skin recovers.
- It is forbidden to apply the device on sensitive parts or private parts (such as eye, lip, nipple, vagina, etc.).
- People with polytrichia or other serious complications are not recommended.
- People with sensitive skin should pay attention to the ice-cool function. Stop using this product if you don't feel comfortable and consult your doctor for advices.
- In winter, people with dry skin should pay attention to daily skin care like maintaining moisture.
- Please wear goggles when operating. Do not directly stare at the flashing light, which would do great damages to the eyes.
- Though slight redness and itching of the skin rarely occurs, it is normal phenomenon, if so, wait for the skin to cool itself.
- Do not spill alcohol directly on the flashing window when cleansing. Please use cotton swabs or dust-free cloth dipped in a small amount of medical alcohol (or use alcohol cotton) to clean the stain.
- Avoid using gel cosmetics or liquids when operating the devices so as not to damage the flashing window.
- Avoid eating food like spinach, rape or mustard during the operation.
- Adjust the gear accordingly; if the skin allows, the higher the gear the better the effect, people with darker skin color should be careful when choosing high gear; otherwise it is easily to get the skin burned.

- After usage, do not touch hot water within 8 hours and pay attention to skin humidification and sunscreen.
- Avoid ultraviolet and computer radiation during usage and wear sunscreen when going out.
- During usage, the edge of flashing window may get stained or black spot because of the residual hair or dirt left on the edge. Although this doesn't affect functions of the device, please keep the flashing window clean at all time.
- This device is not water proof, so please do not use it in moisture environment like shower room.
- If any malfunction takes place, send the device back to manufacturer to confirm whether it can be repaired. Safety of people always comes first under any circumstances.
- Ensure that the device is far from sites with strong magnetic walls, moisture, heavy dust, or high temperature.
- Do not make a strong impact or knock on this device. It is prohibited to use it with purposes other than designated care.
- Besides the authorized personnel, others shall not open up the device; otherwise there will be a risk of electric shock. It is also strictly forbidden for non-professionals to carry out maintenance without authorization.

8. Common problems and solutions

Problems	Reasons	Solutions
Unable to boot up	The power supply is not connected	Make sure the adapter or power cord is plugged in
The adapter is hot in use	Normal phenomenon	No actions required
Stains found on the edge of flashing window	Hair or dirt left on the edge and get burned during usage	Use alcohol cotton to clean up the edge before usage
Effect is not satisfactory	Gear is too low; Usage is not enough	Adjust to suitable gear; Follow the guidance strictly
No illumination	Flashing window is not fully covered by skin	Make sure the window forms an angle of 90 degrees on the skin
Automatic overheating protection(Error 3)	Self-protection program for the device being overheated	Put the device in a ventilation place to cool down
Illustration causes allergy	Gear is too high; Flashing window is damaged	Adjust to suitable gear; Contact customer service for help
Peculiar smell comes with illustration	The flashing window is covered by dust; There is hair left on skin	Use alcohol cotton to clean up the edge; Make sure the hair is shaved to not release any residual left on skin