


MONRO
Multifunctional Anti-Aging Skin Care Device

STRUCTURE COMPONENT

INTRODUCTION

CB12
Instructions for use
Please read the instructions carefully before use

www.praimy.com hello@praimy.com

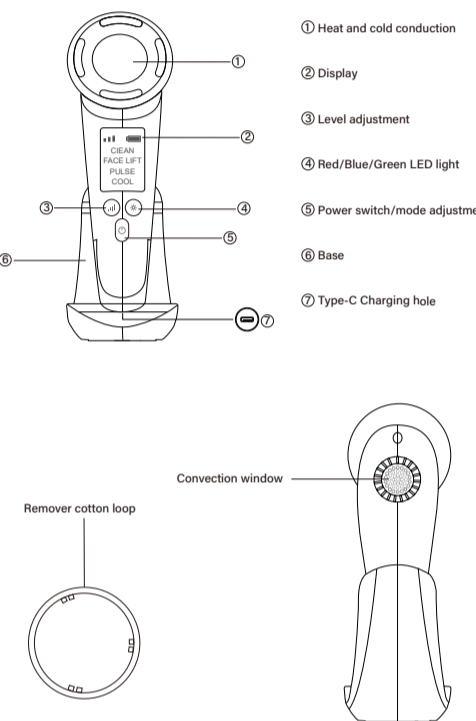

Remover cotton loop

Convection window

1. Micro vibrations beauty principle

Each person's face has at least three million pores, while cleansing liquid and hands wash only clean skin surface, it is difficult to remove the oil pores, only clean out the wash only clean skin surface, it is difficult to remove the oil pores, only clean out the 46 years specialize in skin care experience, the whole development process guided by skin specialists, combined with our advanced engineers to micro-vibration technology, fitting face of intelligent design, all direction clean your skin, 8000 times per minute, high-frequency vibration, vibrate out of deep pores oil and deeply clean every pores. Pores lightweight transparent, delicate skin smooth, once again change the skin more through cleanliness.

2. Principles of hot and cold beauty

4. Japanese professor shows that repeatedly reaching the skin at a temperature of 42 °C helps to produce protein, repair skin cells and strengthen skin elasticity to anti wrinkles. The warm 42°C treatment helps to open the pores for deep cleaning, while the cold treatment with a cold effect of it constantly tightens pores and locks in nutrition. Alternating warm and frozen skin care, it stimulates the liquid circulation at the bottom of the skin and speeds up the metabolism of the cells. Besides effectively reducing the elimination of dark eye circles and simple eye lines, it also brightens skin. Even without makeup, the skin is also be shiny and bright. Combined with micro-vibration, it is beneficial to enhance the skin's firmness and relieve fatigue.

Ice sense effect

- Help skin firmness, maintain flexibility, prevent allergies and other problems
- Shrink pores, lock moisture and anti-wrinkle
- Can soothe muscles, relieve skin fatigue.

Warm effect

- Warning: makes people feel warm and, helps to relax muscles, relieve pain;
- Putting a sponge to deep cleansing, better able to expor the fa cetoxins, blackheads, dirt.

3. EMS beauty principle (Electrical Muscle Stimulation)

EMS is used to activate the skin cell, stimulating rejuvination of collagen, to improve skin elasticity, soften wrinkles and leave skin looking smooth and healthy.

4. Intermittent pulse

The super pulse current is released intermittently, which can achieve the acupuncture feeling of tapping and soothing the muscles. First use during care can make the effect more prominent.

5. LED photon beauty principle

- Red light wavelength 630nm(±10nm) deeply reaches into the 630nm skin to promote blood circulation and the collagen regeneration enhance the skin elasticity, to make the skin to be smooth and moisture. Be suitable for the slack & wrinkly skin.
- Blue light wavelength 470nm(±10nm) deeply reaches into the 480nm skin, to kill the bacteria, promote the blood circulation, can penetrate into 0.25mm under the cortex, destroy the bacteria living environment, anti-inflammatory, bactericidal, remove acne, acne scars and so on effectively.
- Narrow general green light 550nm(±10nm) can repair damaged skin, it can be used during regular massage treatments or in directing and exporting. Green light can brighten the skin's luster and transparency, inhibit the production of melanin and promote the production of collagen, Speed up the skin metabolism cycle.

SPECIFICATION

OPERATION INSTRUCTIONS IN DETAIL

SAFETY INSTRUCTIONS

WARRANTY CARD

Name: Multifunctional Anti-Aging skin Care Device
Model: MONRO
Battery capacity: 1000mAh
Input voltage: DC 5V
Temperature range: -6-42°C(±5°)
Vibration frequency: 8000(±10%)rpm
Charging time: approximately 3H
recharging current: 480MA

Net weight of host: 162g
Net weight of base: 71g
Work time: 5min/Time or 2min/Time
Host size: 176*60*50mm
Use power: 6W

3. Press the power button to switch to Pulse mode. Hold the instrument in contact with the electroplated decorative parts on the side, and at the same time the message head is close to the face, the instrument will vibrate. And there is a finger-pressing sensation of tapping and tapping, according to personal needs, you can press to button ④ adjust temperature and pulse intensity, the intensity is divided into three levels, each mode can be adjusted by pressing button ⑤ adjust three light modes (5 minutes smart timer).

4. Press the power button to switch to cool mode, the beauty device starts to cool, hold the device on the contact side. The surface is electroplated decorative parts, while the message head is close to the face, the instrument will vibrate. According to personal needs you can press to button ④ adjust the temperature range, intensity is divided into three levels, each mode can be adjusted by pressing button ⑤ the red light and blue light mode (Smart timer for 2 minutes, the fan will still work for 30 seconds to dissipate heat after shutting down).

2. Using Method as following pic and text:

1 Clean mode

Need to be used with a cotton pad, attach the cotton pad, pour the makeup remover, turn on the clean mode, warm it up + Positive ion export + vibration massage can effectively clean the pores and clear out (According to the needs of personal skin type ④ you can press the button to adjust the three kinds of lights)

2 EMS Lifting mode

After cleansing, apply skin care products on the face, turn on the EMS lifting mode, EMS micro current stimulates muscles + absorption/massage makes the skin more vigorous and elastic, improves the fine lines on the skin surface, and restores the skin's youthful, smooth and soft soft and white. (According to the needs of personal skin type ④ you can press the button to adjust the three kinds of lights)

1 Pulse mode

After cleansing, apply skin care products on the bottom on the pulse mode, intermittent pulse current + vibration massage can soothe the skin with high intensity, strengthen the introduction effect and direct the pulse to the deep skin to reduce wrinkles. (According to the needs of personal skin type ④ you can press the button to adjust the three kinds of lights)

2 Cool pack mode

Apply the mask on the face, turn on the cool compress mode, cool compress + vibration massage can effectively shrink the pores and lock skin moisture. (According to the needs of personal skin type ④ you can press the button to adjust the three kinds of lights)

Tips

- The power consumption is relatively high at cool pack mode. After shutting down, the cooling fan will remain on for a few minutes to dissipate working heat. Please wait for at least 10 minutes before the next use.
- If the vibration is intermittent during use, it may be due to dry skin or poor contact with the instrument. It can increase the moisture of the skin properly, and use it tightly and close to the skin.
- When the instrument under using, if there is a "beep" sound and the battery symbol flashes, it means that the instrument's power is low and it needs to be recharged.
- The battery symbol flashes when the instrument is charging, and keep the lights on when it is fully charged.

1. Please read the instructions carefully before use, and then use this product correctly in accordance with the operating instructions.

2. Please do not disassemble or modify the machine by yourself. If the product fails, please consult the relevant customer service of the purchase store.

3. Do not pull down or remake the device, please consult the distributor for repairing.

4. Only use the device for its intended purpose as described in this manual.

5. Do not store or use the device in the high temperature and humidity place such as under sunshine or near fire or in the bathroom, etc. To avoid fire or electric shock. Do not use or place it in near the high temperature appliance.

6. Do not use and place the device in the humid condition, in case of any short circuit, electrical shock or fire.

7. Keep away from children and physical disabled people allow to use/lose supervision is necessary when this device is operated by use under 18 years of age.

8. Please don't let children or people with limited physical activity use it. Use under supervision.

9. Please don't let children or people with limited physical activity use it. Use under supervision.

10. After using the ice pack mode, the message head of the machine head will heat up, please do not use it continuously.

WARNING:

- Please stop using when product abnormal heating to avoid burns.
- Please do not drop or clash the device.
- Please apply skin care products during use to avoid dry skin.
- When the Type-C data cable is damaged, please do not use it to avoid short circuit and fire.
- Below conditions are not allowed to use:
Heart disease, fever, communicable disease, cancer, open cuts, acute disease, hypotension, hypertension, install pacemaker, artificial heart or lung, during medical

1. The free warranty of this warranty card is limited to electrical performance failures caused by product, technical and quality problems.

2. One year guarantee from the date of purchase.

3. Free to change broken parts, freight is collected by customers.

4. If breakdown caused by the human factor, free maintain won't be given and freight shall be collected by customers.

5. Please fill the guarantee clearly, without any revision. Please take good care of guarantee. If not, won't be re-provided.

Product Name	Item No		
Customer Name	Phone		
Purchase Date	Means Of Purchase		
Address			
Records			
Date	Content	Results	Service Provider

Implementation of standards: GB 4706.1-2005
GB 4706.1-2008

