

Thank you for purchasing this product! In order to use this product correctly, please read this manual before use; in order to read the manual at any time when you need it, please keep this manual in a safe place.

About this manual:
The product pictures and related product features shown in this manual, as well as some nouns, can only be used as a reference. The user will use the actual product as the standard, and the referenced data comes from the company's laboratory.

Safety Precautions

Be sure to read the safety precautions before use. Please use this product correctly according to the contents and recommendations of the manual. The product should be cleaned in time after use to avoid accumulation of grease and its residues.

	Related content: Incorrect actions can lead to personal injury
	Related content: Incorrect operation can result in personal safety or property damage

- Danger**
 - Please use the charging cable supplied with this product or the charging cable that meets the product specifications, and connect the charger plug to the household socket correctly to avoid malfunction or fire caused by heat.
- Since this product has a built-in rechargeable battery, do not put it in a fire or use it in a high temperature environment, which may cause malfunction or fire due to heat.
- If the product is unusually hot or deformed, damaged, or smelly, do not use it to avoid the risk of electric shock and fire.
- It is strictly forbidden to use a person who has undergone cardiac surgery or a medical electronic device such as a pacemaker or an embedded defibrillator placed inside the body.




- To prevent the occurrence of danger, do not use it with medical electronic equipment.
- For first-time use or sensitive skin, please try a small area of skin and confirm the condition of the skin before using.

1

- If the machine is damaged, in order to avoid danger, it must be repaired by the manufacturer, its maintenance department or a specialist in a similar department.
- Do not use this product (such as alcohol, benzene, spray, etc.) near flammable or explosive materials, which may result in fire or explosion.
- Place packaging materials such as plastic bags where children can't reach them.
- Please turn off the power before cleaning and maintenance.
- Do not use this product for purposes other than those.
- When using, if the product is abnormally hot, please stop using it immediately and shut it down.
- Do not use the machine if it is faulty. Using this product in a faulty state may have adverse effects on the body or other hazards.
- Please do not use this product continuously in the same area, otherwise it may cause a burning sensation.
- This product is not intended for use by persons with physical, sensory or mental problems, or those who lack experience and knowledge (including children), unless supervised by a guardian and instructed in the use of the product, and responsible for the safety of the user.
- Do not immerse the appliance in water or rinse it directly with a water tap.
- If you feel any discomfort in your skin or body, stop using it immediately.
- Please do not store this product in a humid environment.
- When using, remove precious metal objects and watches.
- Do not use this product in the following situations:
 - During pregnancy, lactating or menstruating women.
 - Allergic people have skin diseases, skin inflammation, skin irritation, other skin diseases and severe sunburn.
 - People with skin damage, wetness, swelling, etc.
 - The person who has just finished surgery and is recovering.
 - Due to the upper fracture or the plastic face, the metal has been received on the face or in the body, and the silicon polymer is embedded in the surgery.
 - People with heart disease, malignant tumors, epilepsy, diabetes, infectious diseases, acute conditions.
 - Long-term use of steroid hormone drugs or liver dysfunction caused by capillary vasodilatation.
 - In the past 3 months, surgery has been performed at the site of use, with long-lasting laser or deep chemical stripping.

2

Part Identification

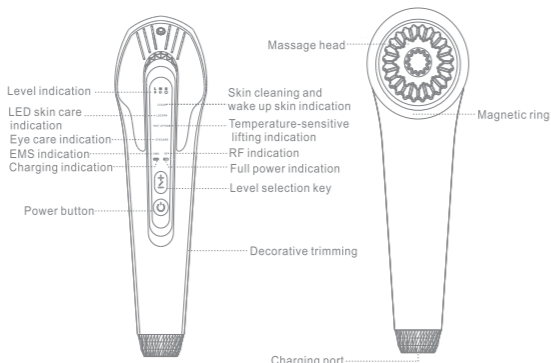

【Description on Panel】

- Power on: short press the power button; power off: long press the power button for 2 seconds.
- 10 minutes automatic shutdown function.
- Mode switching: The default is "cleaning" mode, each time you press the power button to switch one mode, a total of 6 modes (cleaning mode, LED skin rejuvenation mode, warming lifting mode, eye care mode, EMS Mode, RF mode.)
- Gear selection: The default is 1 gear at start-up, each time you press the gear selection key, it will change 1 gear, a total of 3 gears.


Operation

【Charging】

For the first time, please fully charge the machine before using it.

3

- Connect the machine to the charging cable and charger, and plug in the household power socket;
- During charging, the charging indicator flashes;
- When fully charged, the battery indicator is always on.
- Charging time takes about 3 hours;
- Fully charged up to 3 times (10 minutes each time);
- The number of uses will vary somewhat depending on usage.

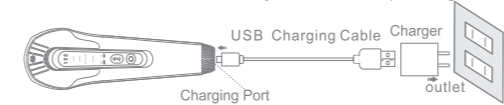

- ※ Do not use on Adam's apple, eye balls, mucosa and more positions.
- ※ During using, do not stay in a certain position or press the skin hardly.
- ※ For the first time use, it is advised to start from the 1 level, and then choose proper strength levels according to your needs.
- ※ When using, the hand must touch the decorative trimming

【Preparation before use】

- Prepare skin care products such as essences, masks or creams.
- When taking care around the eyes, remove the contact lenses first.
- Remove makeup according to the daily makeup removal procedure.
- Clean the massage head with a moist cotton pad.

【CLEAN】 Skin cleaning and wake up skin

- Remove magnetic ring with fingers.
- Put cotton pad on head and fasten it with magnetic ring
- Pour toner on the cotton pad
- Hold the main machine by hand and carry out facial care with corresponding caring method.
- After the care is finished, press the Start button for 1 second to turn off the apparatus.

4

- Press power button, defaults is CLEAN mode
- Press the Level button to select the appropriate level.
- Moves slowly from the inside out, from the bottom up to clean face

- After the care is finished, press the Start button for 1 second to turn off the apparatus.
- Clean face with clean water after cleansing
- Remove the dirt cotton pad and wipe the massage head clean.

※ Note: clean mode is recommended for 2 to 3 days.

【Facial Care】

- Apply a mask or evenly spread daily use lotion on your face.
- After starting, press the Start button to select the relevant mode.
- Press the Level button to select the appropriate level.

- Hold the main machine by hand and carry out facial care with corresponding caring method.
- After the care is finished, press the Start button for 1 second to turn off the apparatus.

5

- Follow the contour of the face do upward pulling action.
- From the inside of the ear to the back of the neck, from top to the bottom.

- From the eyebrows to the left and right sides to do a symmetrical pull action.
- Pulling from the forehead to the hairline.

- Move around the outer corner of the eye and move up from the bottom of the eye to the outer corner of the eye.
- From the outer corner of the eye, consciously lift the skin for maintenance.

- Starting from the nasolabial folds, along the expression muscles from the bottom up to the corners to do the lifting action.
- Doing a circular motion

- Lifting upwards from the nasolabial folds to cheeks
- Lifting upwards from chin to cheeks

6

Tips

- During the nursing process, the skin may feel hot or slight pain. Some users may feel a slight redness after using the skin. This phenomenon usually disappears within a few minutes after use. If the skin feels abnormally hot during use, it is recommended to stop using and let the skin cool down first.

Maintenance

- Protect the appliance from moisture and do not place the appliance in a dishwasher or microwave oven.
- This product is not waterproof and should not be immersed in water.
- Do not use sharp or corrosive, volatile thinner cleaners.
- It can be cleaned with a soft damp cloth and wiped with a soft, dry cloth.
- Do not use for a long time, please take out the machine every 3 months or so before fully storing it.

Fault analysis and elimination

In case of any malfunction during the operation, please carry out following inspections. If failures cannot be eliminated or the product cannot operate properly, please contact after-sales service center. Do not disassemble and attempt to repair it by yourself.

Situation	Possible Reason	Measure
Sudden shutdown or machine does not start	Low battery	Fully charge the host and use it again
Charging light is not lit	The main unit, charging stand, and charger are not connected	Reconnect the main unit and charger, or plug the charger into another household socket.
	Battery is dormant	The battery needs 6-30 minutes back to working. It is a normal phenomenon and needs to be fully charged before use.
Massage head contact with skin without irritation or feeling of heat.	The skin is too dry	Please apply the lotion to your skin before use

7


After-sales

- This appliance is intended for personal use only and is not suitable for commercial use.
- We offer a one-year warranty on the product, which is calculated from the date of purchase. During the warranty period, we will replace parts or the entire product due to the material or the production process.
- To apply for a warranty, the product must be shipped to an authorized repair station or factory along with a
 - Damage caused by improper use, such as using the wrong power supply, inserting the product into an improper socket, man-made damage, etc.
 - The product caused by normal use is aging and wear, but it has little effect on the function of the product.
 - Product problems caused by other uses, or improper handling or non-professional repairs or the use of accessories other than the company.
 - Consumables or bonus accessories.
 - Damage caused by force majeure factors (such as lightning, earthquakes, fires, floods, etc.).
- After-sales service does not include damage caused by the above reasons.

Waste Treatment


The marking on this product or packaging material or accessory indicates that the material is made from a recyclable material. Separate waste into different recycling levels to facilitate the recycling of useful materials. The lithium battery is installed inside the appliance, so:

- before removing the utensils, the battery must be removed from the appliance.
- when the battery is removed, the appliance must be shut down for power failure.
- batteries should be disposed of safely.
- Users cannot replace batteries by themselves.

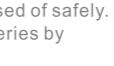

8

- Steps to remove the battery:
- Open the decorative cover and unscrew the screws above.
 - Pry the product rod cover with a tool and remove it.
 - Unplug the battery leads and remove the battery.

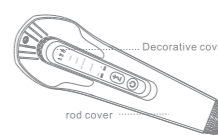

Specification

Product name: Microcurrent Face Lifting Device
 Model: CLEOPATRA
 Rated Power: 10.5W
 Charging Power Supply: 5V 1A
 Product Dimensions: 185x55x45mm
 Packing Dimensions: 216x95x60mm
 Packing List: 1 main unit, 1 charging cable, 1 operating manual.

Note: In view of the consistent policy of continuous improvement of products, the company reserves the right to change the product specifications, the actual goods are subject to shipment.

9

CLEOPATRA Microcurrent Face Lifting Device



www.primy.com

hello@primy.com