

Thank you for purchasing this product! In order to use this product correctly, please read this manual before use; in order to read the manual at any time when you need it, please keep this manual in a safe place.

About this manual:

The product pictures and related product features shown in this manual, as well as some nouns, can only be used as a reference. The user will use the actual product as the standard, and the referenced data comes from the company's laboratory.

Safety Precautions

Be sure to read the safety precautions before use. Please use this product correctly according to the contents and recommendations of the manual. The product should be cleaned in time after use to avoid accumulation of grease and its residues.

	Related content: Incorrect actions can lead to personal injury
	Related content: Incorrect operation can result in personal safety or property damage

Danger


- Please use the charging cable supplied with this product or the charging cable that meets the product specifications, and connect the charger plug to the household socket correctly to avoid malfunction or fire caused by heat.
- Since this product has a built-in rechargeable battery, do not put it in a fire or use it in a high temperature environment, which may cause malfunction or fire due to heat.
- If the product is unusually hot or deformed, damaged, or smelly, do not use it to avoid the risk of electric shock and fire.

Caution

- To prevent the occurrence of danger, do not use it with medical electronic equipment.
- It is strictly forbidden to use a person who has undergone cardiac surgery or a medical electronic device such as a pacemaker or an embedded defibrillator placed inside the body.
- For first-time use or sensitive skin, please try a small area of skin and confirm the condition of the skin before using.

- If the machine is damaged, in order to avoid danger, it must be repaired by the manufacturer, its maintenance department or a specialist in a similar department.
- Do not use this product (such as alcohol, benzene, spray, etc.) near flammable or explosive materials, which may result in fire or explosion.
- Place packaging materials such as plastic bags where children can't reach them.
- Please turn off the power before cleaning and maintenance.
- Do not use this product for purposes other than those specified in this manual to avoid malfunction or fire due to heat.
- When using, if the product is abnormally hot, please stop using it immediately and shut it down.
- Do not use the machine if it is faulty. Using this product in a faulty state may have adverse effects on the body or other hazards.
- Please do not use this product continuously in the same area, otherwise it may cause a burning sensation.
- This product is not intended for use by persons with physical, sensory or mental problems, or those who lack experience and knowledge (including children), unless supervised by a guardian and instructed in the use of the product, and responsible for the safety of the user.
- Do not immerse the appliance in water or rinse it directly with a water tap.
- If you feel any discomfort in your skin or body, stop using it immediately.
- Please do not store this product in a humid environment.
- When using, remove precious metal objects and watches.
- Do not use this product in the following situations:
 - During pregnancy, lactating or menstruating women.
 - Allergic people have skin diseases, skin inflammation, skin irritation, other skin diseases and severe sunburn.
 - People with skin damage, wetness, swelling, etc.
 - The person who has just finished surgery and is recovering.
 - Due to the upper fracture or the plastic face, the metal has been received on the face or in the body, and the silicon polymer is embedded in the surgery.
 - People with heart disease, malignant tumors, epilepsy, diabetes, infectious diseases, acute conditions.
 - Long-term use of steroid hormone drugs or liver dysfunction caused by capillary vasodilatation.
 - In the past 3 months, surgery has been performed at the site of use, with long-lasting laser or deep chemical stripping.

Part Identification

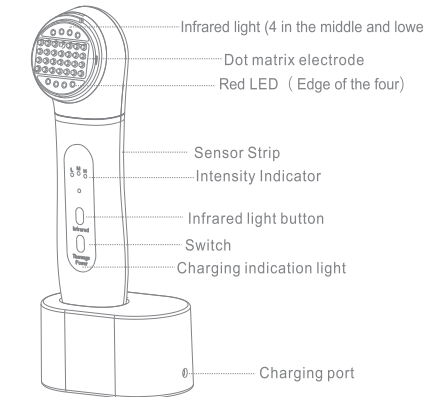

How to Use


[Charging]

For the first time, please fully charge the machine before using it.

- Connect the machine to the charging cable, charger and plug in the household socket;
- During charging, the charging indicator lights up in red;
- When fully charged, the charging indicator lights up in green;
- Charging time takes approximately 3 hours.
- Fully charged for 3 hours;
- Depending on usage, the usage time will vary.

[Preparation before Use]



- Wipe the machine clean with a soft cloth or cotton pad;
- Please remove the metal objects (rings, necklaces, earrings, bracelets, etc.) from your body.
- Please remove makeup and cleanse before use to ensure that your skin is in good condition.

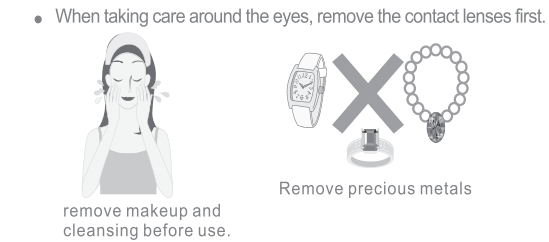

Tips:


- This product should be used together with nutrient solution or beauty gel.
- Please do not use this product on the throat, upper eyeball, mucous membranes, etc.
- Do not stop in one place or squeeze the skin while using it.
- For the first time, it is recommended to start from the L intensity, and gradually select the appropriate gear according to your personal skin feeling.
- When used, the sense of current is based on the skin's own moisture level. After a set of processes, the skin is hot. When the care was first started, the skin did not feel the obvious heat. After a few minutes, it began to slowly heat from the inside out.
- During the nursing process, the skin may feel hot or slightly tingling. Some users may feel a slight redness after using the skin. This phenomenon usually disappears within a few minutes after use. If the user feels abnormal skin fever during use, it is recommended to stop using it and first cool the skin.
- For the first half month of the Thermage using, it is recommended to use it once every other day and then every 4 days.

[Operation]

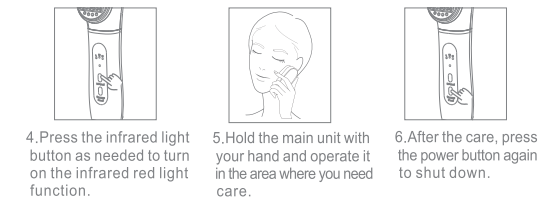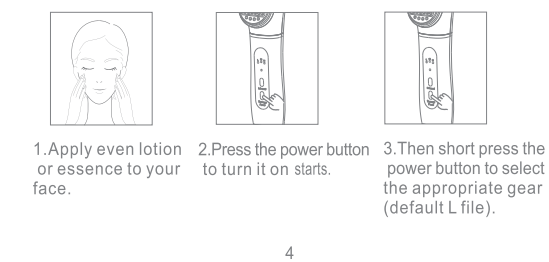


[Nursing Method]



Repair key parts such as wrinkles, marks and spots

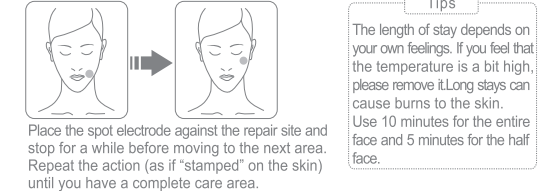













Leading-in

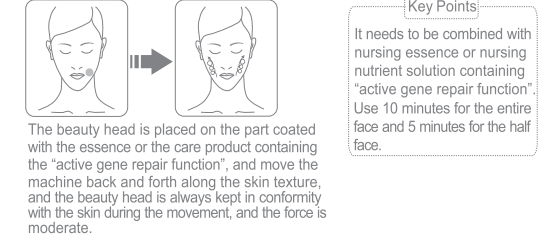

Dredge lymphatic care

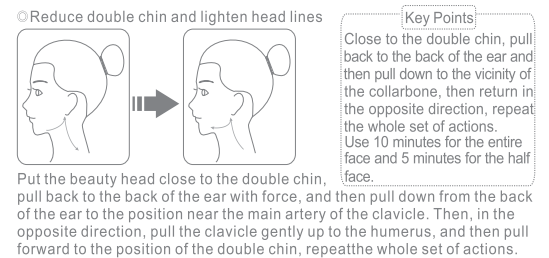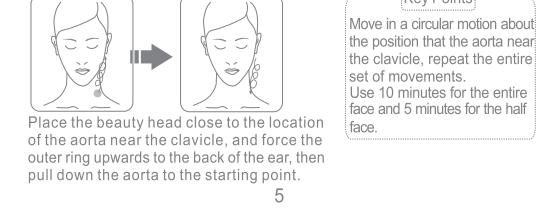

Reinventing the beautiful V-shaped face

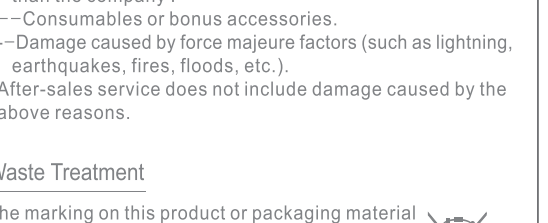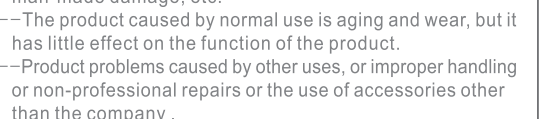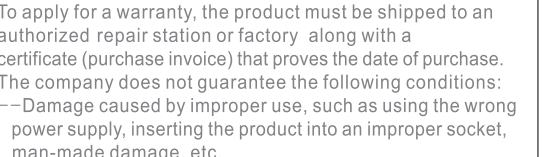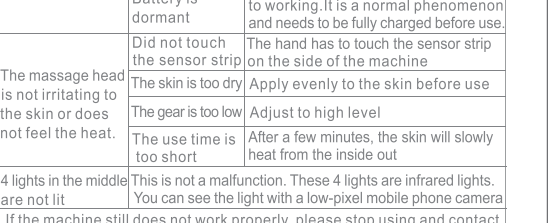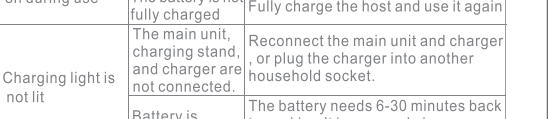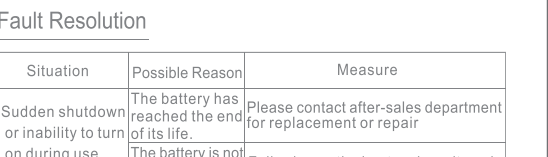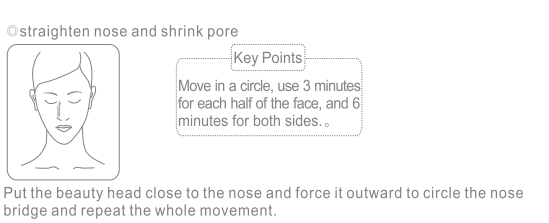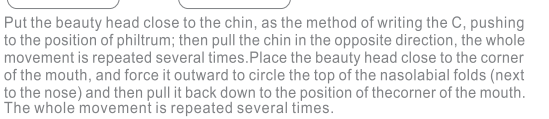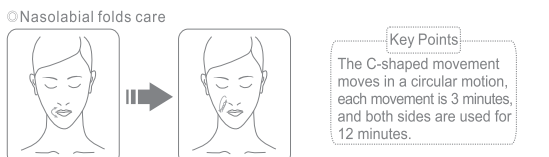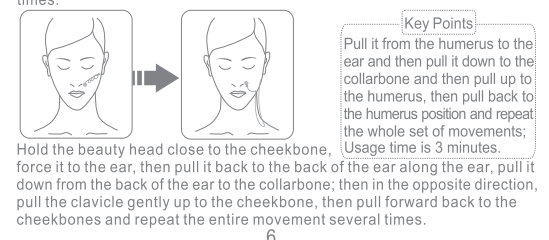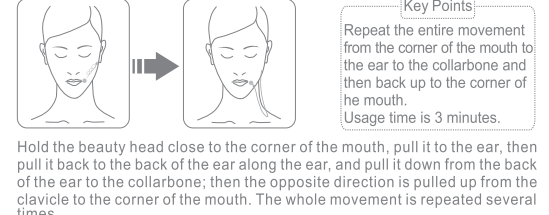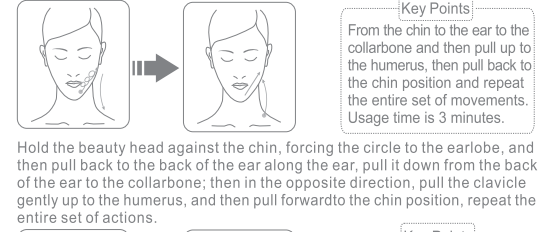

Fault Resolution

Situation	Possible Reason	Measure
Sudden shutdown or inability to turn on during use	The battery has reached the end of its life. The battery is not fully charged.	Please contact after-sales department for replacement or repair. Fully charge the host and use it again.
Charging light is not lit	The main unit, charging stand, and charger are not connected. Battery is dormant.	Reconnect the main unit and charger or plug the charger into another household socket. The battery needs 6-30 minutes back to working. It is a normal phenomenon and needs to be fully charged before use.
The massage head is not irritating to the skin or does not feel the heat.	Did not touch the sensor strip. The skin is too dry. The gear is too low. The use time is too short.	Adjust to high level. Apply evenly to the skin before use. After a few minutes, the skin will slowly heat from the inside out.
4 lights in the middle are not lit	This is not a malfunction. These 4 lights are infrared lights, you can see the light with a low-pixel mobile phone camera. If the machine still does not work properly, please stop using and contact the after-sales department.	

Maintenance

- Protect the appliance from moisture and do not place the appliance in a dishwasher or microwave oven.
- Please turn off the power before cleaning.
- This product is not waterproof and should not be immersed in water.
- Do not use sharp or corrosive, volatile thinner cleaners.
- It can be cleaned with a soft damp cloth and wiped with a soft, dry cloth.

After-sales

- This appliance is intended for personal use only and is not suitable for commercial use.
- We offer a one-year warranty on the product, which is calculated from the date of purchase. During the warranty period, we will replace parts or the entire product due to the material or the production process.

Environmental Protection

The symbol on this product or in the instruction manual or on the packaging indicates that this product cannot be disposed of as general household waste at the end of its useful life. Please consult your local government for recycling according to the classification of waste.


Specification

Product Name: 3 in 1 Combination Beauty Device
Model: ELEANOR
Voltage: 3.7V ---
Rated Power: 4.9W
Charging Power Supply: 5V ---1A
Net Weight: 155g
Product Dimensions: 166x46x46mm
Packing List: 1 main unit, 1 charging stand, 1 charging cable, 1 manual.

Note: In view of the consistent policy of continuous improvement of products, the company reserves the right to change the product specifications, the actual goods are subject to shipment.



ELEANOR
3 in 1 Combination Beauty Device

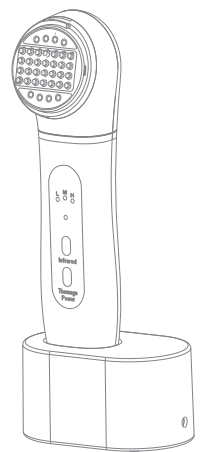

www.pramy.com

hello@pramy.com